

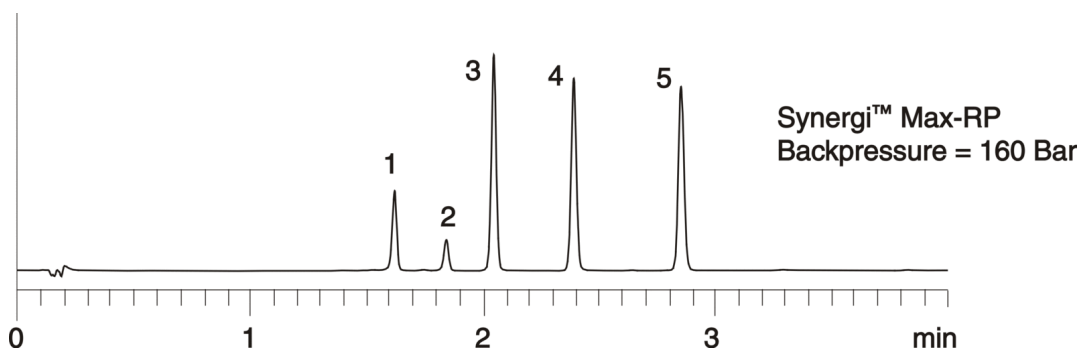
Steroids**Column:** SynergiTM 4 μ m Max-RP 80 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4337-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% Formic acid**Eluent B:** Acetonitrile with 0.1% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	5	95

Flow Rate: 4 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Dexamethasone
- 2 Hydrocortisone
- 3 11- α -Hydroxyprogesterone
- 4 17- α -Hydroxyprogesterone
- 5 Progesterone

