

## Tricyclic Antidepressants on Luna CN

**Column:** Luna® 5 µm CN 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4255-E0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** 20 mM Potassium phosphate pH 4.0

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 0        | 35         | 65    |       |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, CN 4 x 3.0mm, 10/Pk Part No.: AJ0-4305

- Holder Part No.: KJ0-4282



Products used in this application:

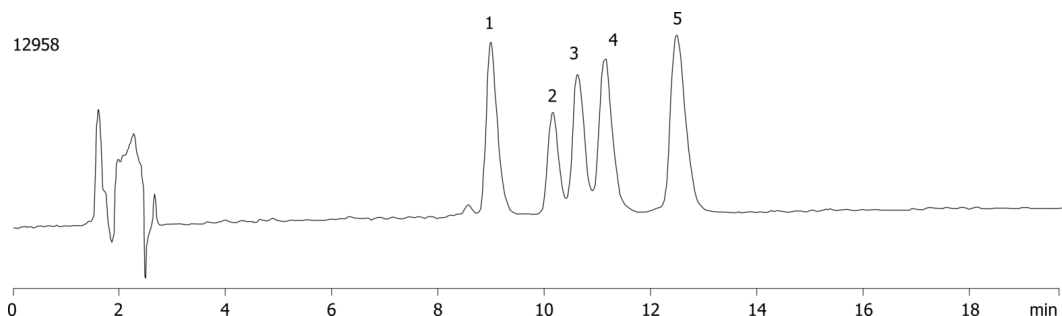

### ANALYTES:

- 1 Doxepin
- 2 Desipramine
- 3 Nortriptyline
- 4 Amitriptyline
- 5 Protriptyline

