

## Isoetharine

**Column:** Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-3014-E0

**Elution Type:** Isocratic


**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (58:35:7)

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 0.7 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 282 nm (ambient)

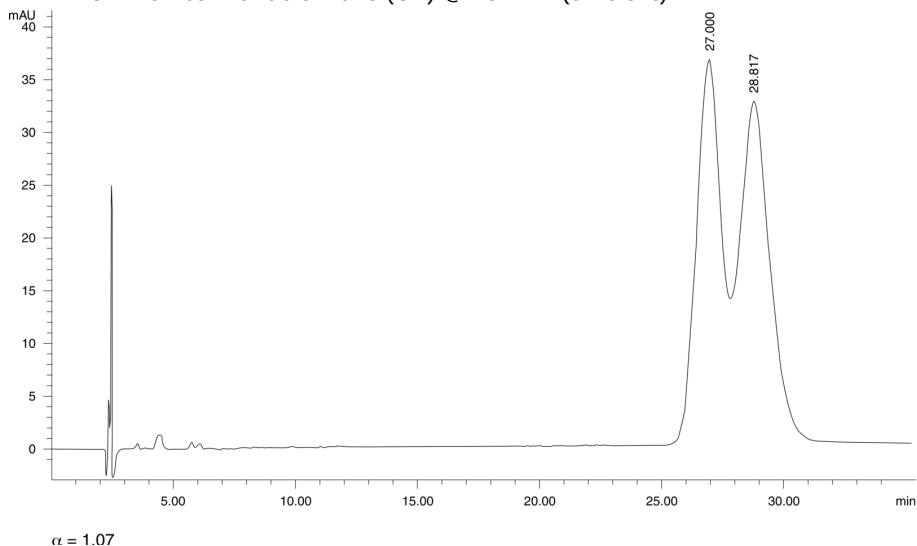

Products used in this application:

**CHIREX®**

### ANALYTES:

**1** Isoetharine

Retention Time: 27 min

**2** Isoetharine

Retention Time: 28.817 min

