


N-Acetylalanine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-3126-E0
Elution Type: Isocratic
Eluent A: 2 mM Copper (II) sulfate in water
Eluent B: Acetonitrile
Gradient


Step No.	Time (min)	Pct A	Pct B
1	0	95	5

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

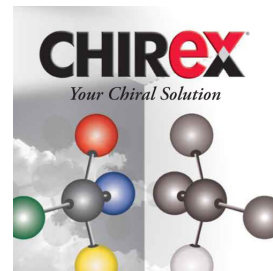

Products used in this application:

CHIREX

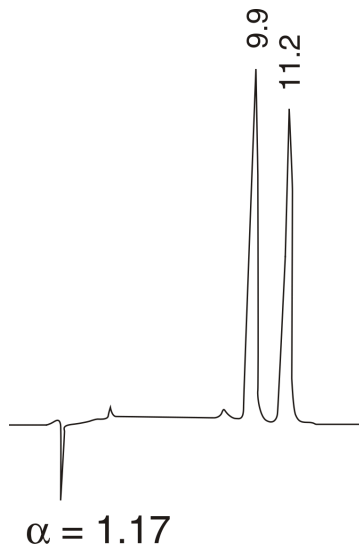

ANALYTES:

- 1 N-Acetylalanine
Retention Time: 9.85 min
- 2 N-Acetylalanine
Retention Time: 11.21 min

