

Normetanephine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 5 mM Copper (II) acetate in water [pH 6.0]

Gradient


Step No. **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 0.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

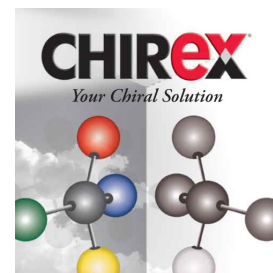

Products used in this application:

CHIREX

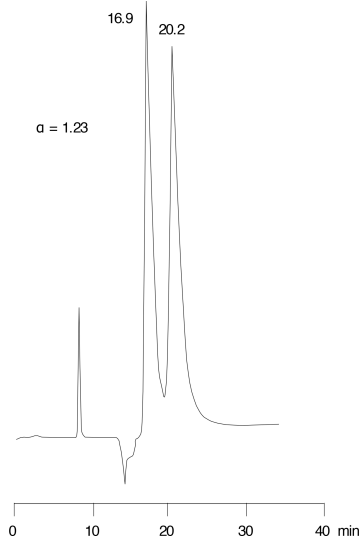

ANALYTES:

1 Normetanephine

Retention Time: 16.9 min

2 Normetanephine

Retention Time: 20.2 min

