

## Basic peptides on Jupiter C18 in 10mM TFA

**Column:** Jupiter® 5 µm C18 300 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4053-E0

**Elution Type:** Gradient

**Eluent A:** Water with 10mM TFA

**Eluent B:** Acetonitrile with 10mM TFA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 90    | 10    |
|                   | 2        | 20         | 55    | 45    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 35 °C

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4321

- Holder Part No.: KJ0-4282



Products used in this application:

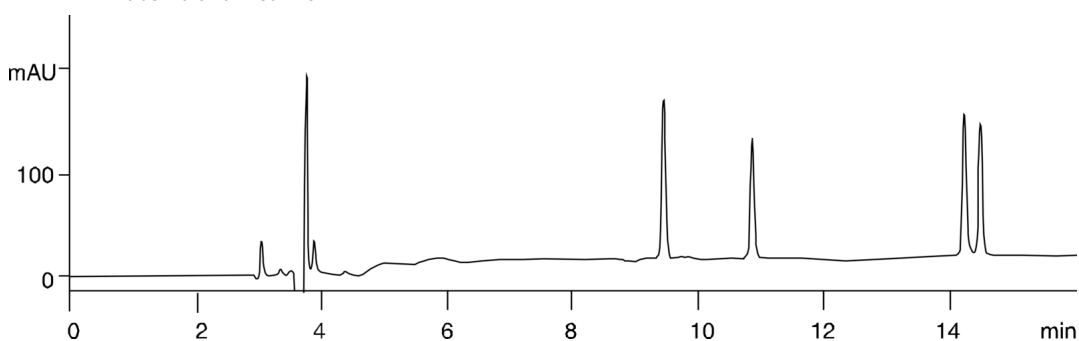

### ANALYTES:

- 1 Basic peptide
- 2 Basic peptide
- 3 Basic peptide 3
- 4 Basic peptide 4

