

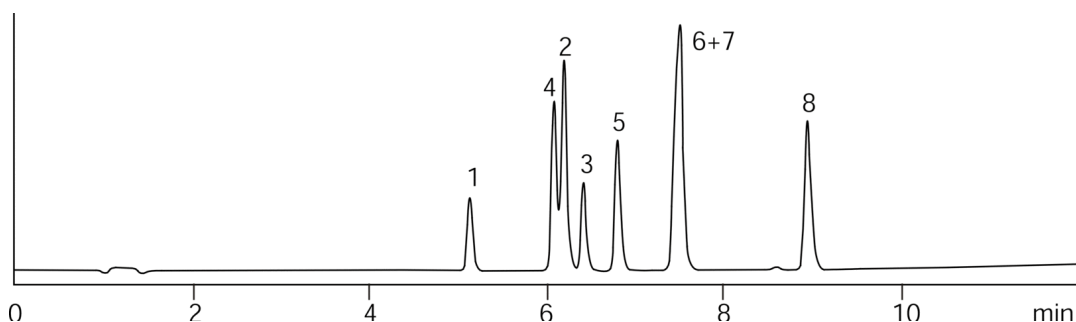
Sedatives on Synergi Max-RP**Column:** SynergiTM 4 μ m Max-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4337-E0**Elution Type:** Gradient**Eluent A:** 20mM Potassium dihydrogen phosphate, pH 4.6**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	15	25	75

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 impurity
- 2 Oxazepam
- 3 Lorazepam
- 4 Clonazepam
- 5 Temazepam
- 6 Chlordiazepoxide
- 7 Desmethyldiazepam
- 8 Diazepam

