

## Phthalates on Luna Amino

**Column:** Luna® 5 µm NH<sub>2</sub> 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4378-E0

**Elution Type:** Isocratic

**Eluent A:** Hexane

**Eluent B:** Ethanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 1        | 0          | 95    | 5     |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, NH<sub>2</sub> 4 x 3.0mm, 10/Pk Part No.: AJ0-4302

- Holder Part No.: KJ0-4282



Products used in this application:

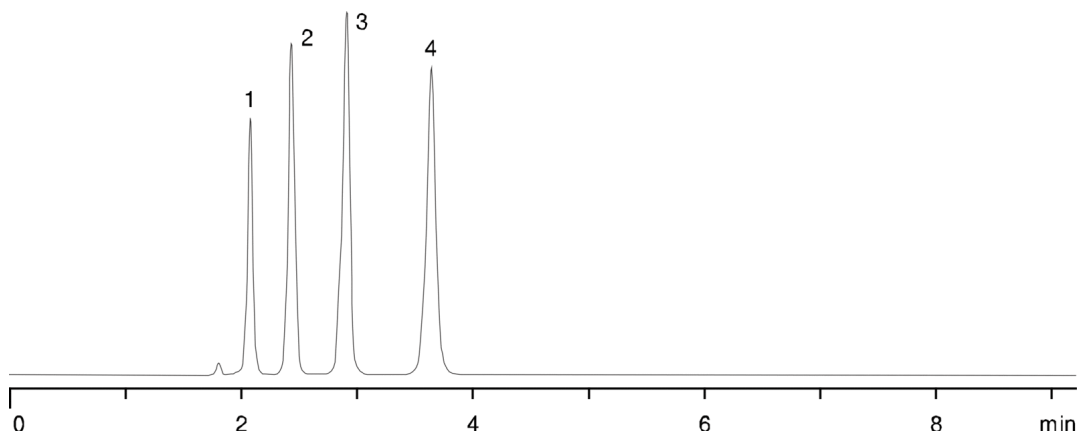

### ANALYTES:

- 1 Diisooctyl phthalate
- 2 Dibutyl phthalate
- 3 Diethyl phthalate
- 4 Dimethyl phthalate

