

## Extraction of Procainamide on Strata™ X

**Column:** Luna® 5 µm C8(2) 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4249-E0

**Elution Type:** Isocratic

**Eluent A:** 20 mM KH<sub>2</sub>PO<sub>4</sub> (pH 7.0)

**Eluent B:** Methanol

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
| 1                 | 0        | 90         | 10    |       |

**Flow Rate:** 1 mL/min

**Col. Temp.:** 25 °C

**Detection:** UV-Vis Abs.-Fixed Wave. (UV) @ 220 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282



Products used in this application:

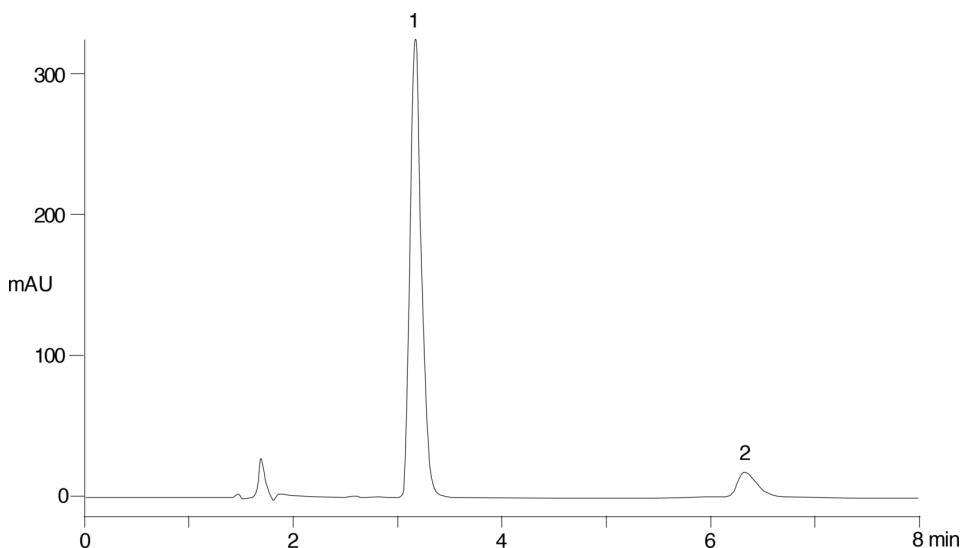

### ANALYTES:

- 1 Sulfanilimide
- 2 Procainamide

