


Abused Inhalants ZB-WAX 30x0.32x0.5**Column:** Zebtron™ ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 µm**Order No:** 7HM-G007-17**Oven Profile:** 40°C for 4 min to 120°C at 10°C/min.**Carrier Gas:** Constant Flow Helium, 3.5 mL/min**Injection:** Split 50:1 100 µL @ 200°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Headspace analysis**Products used in this application:**

Abused Inhalants ZB-WAX 30x0.32x0.5

ANALYTES:

- 1** Hexane
- 2** Ether
- 3** Methyl tert-butyl ether
- 4** Acetone
- 5** Ethyl acetate
- 6** Methanol
- 7** MEK
- 8** Isopropanol
- 9** Ethanol
- 10** Methylene chloride
- 11** Benzene
- 12** Trichloroethylene
- 13** MIBK
- 14** Chloroform
- 15** Toluene

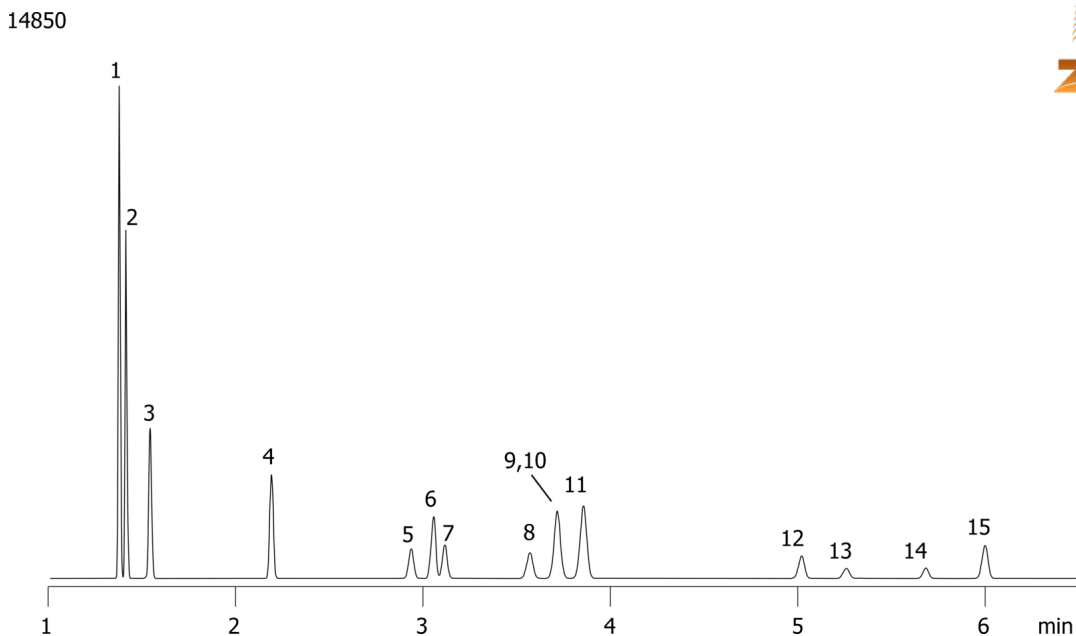