


Nitrosamines ZB-35 30x0.25x0.25**Column:** Zebron[™] ZB-35, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G003-11**Oven Profile:** 40°C to 300°C at 20°C/min**Carrier Gas:** Constant Flow Helium, 0.8 mL/min**Injection:** split 50:1 1 μ L @ 200°C**Detection:** Mass Selective (MSD) (240°C)

14853

**Products used in this application:****ANALYTES:**

- 1 N-Nitrosodimethylamine
- 2 N-Nitrosomethylethylamine
- 3 N-Nitrosodiethylamine
- 4 N-Nitrosodi-n-propylamine
- 5 N-Nitrosopyrrolidine
- 6 N-Nitrosomorpholine
- 7 N-Nitrosopiperidine
- 8 N-Nitrosodi-n-butylamine
- 9 N-Nitrosodiphenylamine

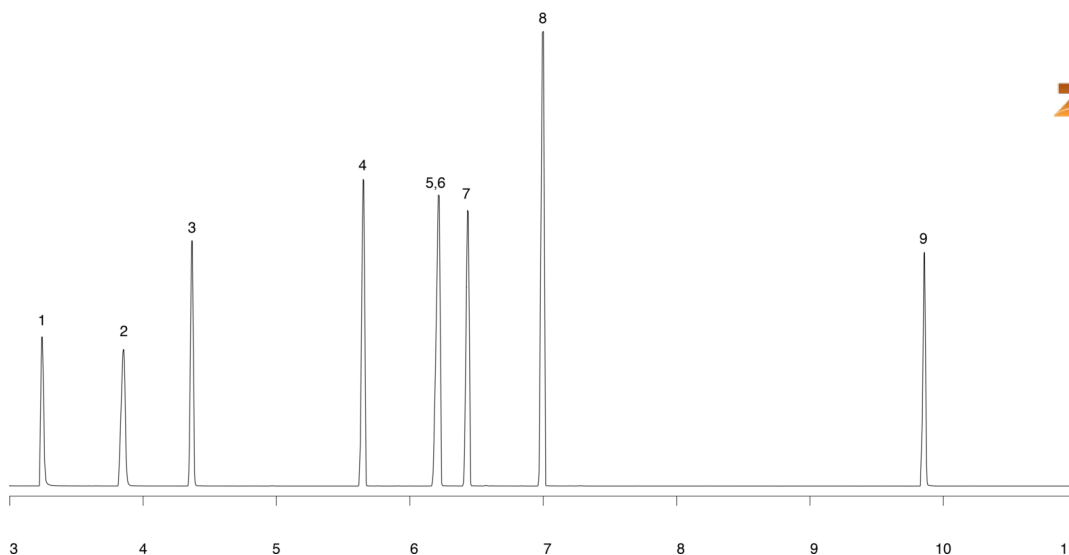