


Aniline - Phenol Separation - ZB-Wax 60x0.32x0.5**Column:** ZebtronTM ZB-WAX, GC Cap. Column 60 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 60 meters x 0.32 mm x 0.5 μ m**Order No:** 7KM-G007-17**Oven Profile:** 130°C to 250°C at 5°C/min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** split 20:1 0.2 μ L @ 275°C**Detection:** Mass Selective (MSD) (240°C)

14861

**Products used in this application:****ANALYTES:**

- 1** Aniline
- 2** Phenol

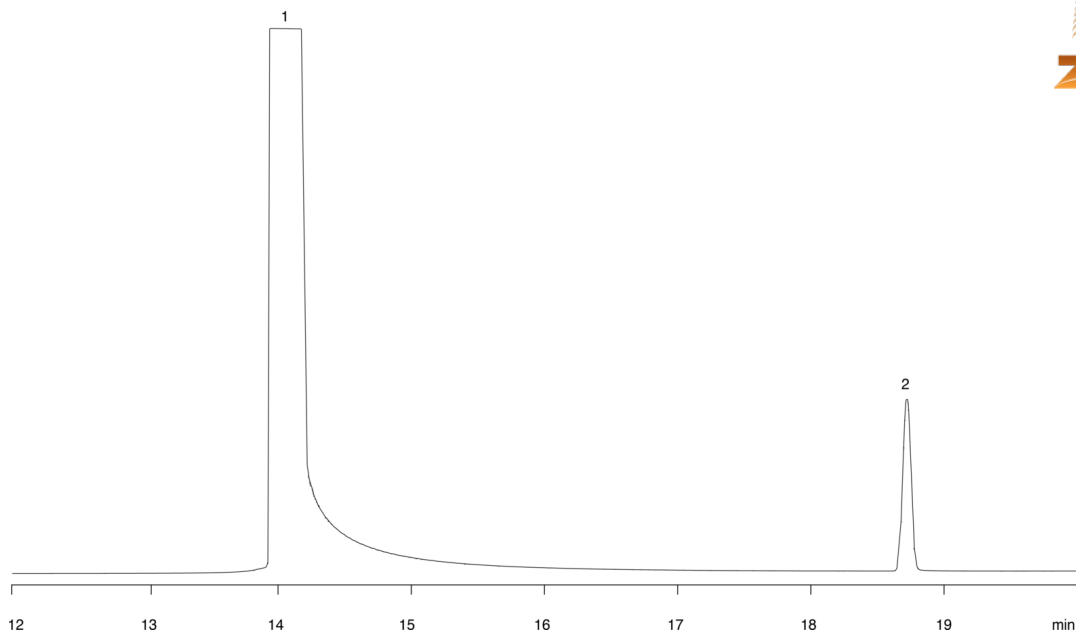