


Phthalate Esters: EPA Method 606 - ZB-5 30x0.53x1.5**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.53 mm x 1.50 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 1.5 μ m**Order No:** 7HK-G002-28**Oven Profile:** 40°C for 6 min to 300°C at 10°C/min for 15 min**Carrier Gas:** Constant Flow Helium, 12.9 mL/min**Injection:** Split 20:1 1 μ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)

14876

**Products used in this application:****ANALYTES:**

- 1** Dimethyl phthalate
- 2** Diethyl phthalate
- 3** Di-n-butyl phthalate
- 4** Butyl benzyl phthalate
- 5** bis(2-Ethylhexyl)phthalate
- 6** Di-n-octyl phthalate

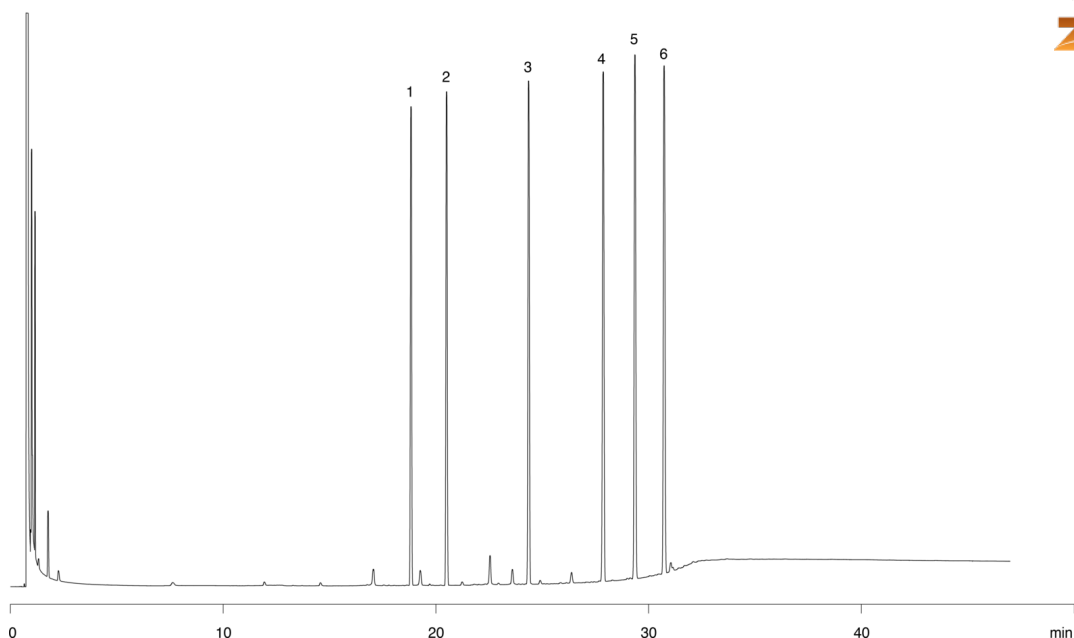