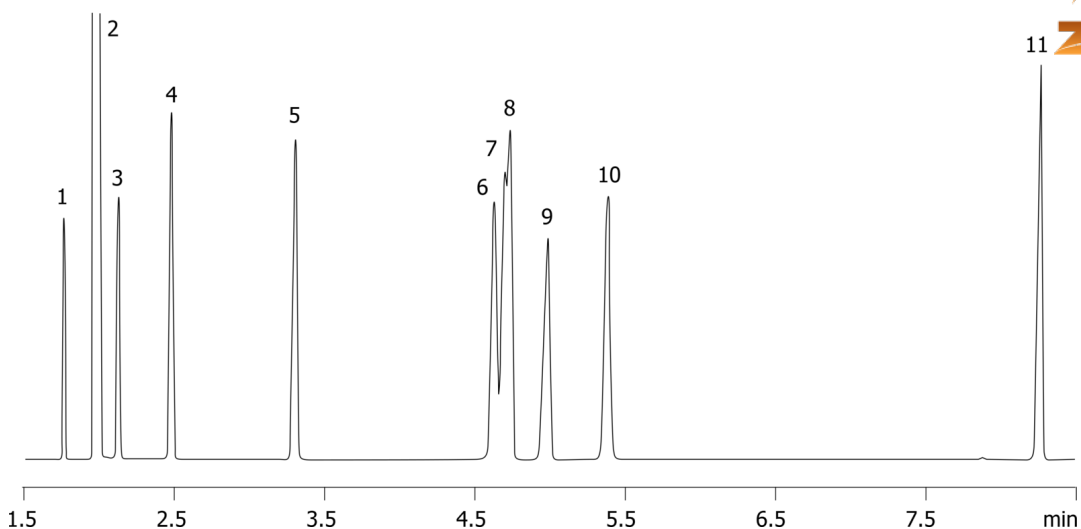


BTEX ZB-50 30x0.32x0.5**Column:** Zebron[™] ZB-50, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** 50% Phenyl 50% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G004-17**Oven Profile:** 60°C to 75°C at 15°C/min to 90°C at 3°C/min for 3 min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** split 20:1 0.2 μ L @ 225°C**Detection:** Flame Ionization (FID) (300°C)

14903

**Products used in this application:**

BTEX ZB-50 30x0.32x0.5

ANALYTES:

- 1 Pentane
- 2 Methylene chloride
- 3 Heptane
- 4 Benzene
- 5 Toluene
- 6 Ethylbenzene
- 7 p-Xylene
- 8 m-Xylene
- 9 o-Xylene
- 10 Decane
- 11 Dodecane

