

IgG Dog Native Acetonitrile

Column: Jupiter® 5 µm C4 300 Å, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4167-B0

Elution Type: Gradient

Eluent A: 0.1% TFA in 95:5 Water/ACN

Eluent B: 0.085% TFA in 95:5 ACN/Water


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	20	5	95

Flow Rate: 0.25 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)

Analyst Note: These applications demonstrate the utility of using low amounts of isopropanol in the organic mobile phase to improve peak shape of hydrophobic proteins. App ID#14918 shows the elution of IgG on a Jupiter 300 C4 column where Acetonitrile is used as the organic mobile phase. App 14918

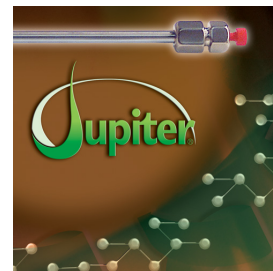

Products used in this application:

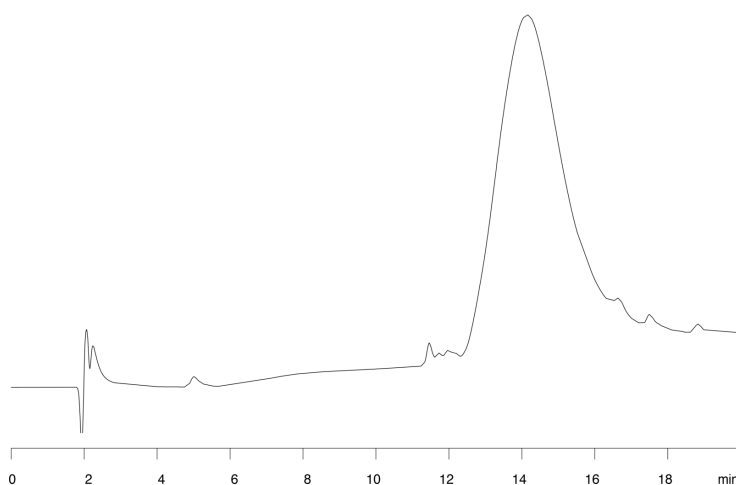

ANALYTES:

1 IgG

Retention Time: 14.14 min

