

## Jupiter C4 P2\_1

**Column:** Jupiter® 5 µm C4 300 Å, LC Column 250 x 2 mm, Ea

**Dimensions:** 250 x 2 mm ID

**Order No:** 00G-4167-B0

**Elution Type:** Gradient

**Eluent A:** 0.1%TFA in 95:5 Water/ACN

**Eluent B:** 60:40 ACN/MeOH

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 80    | 20    |
|                   | 2        | 15         | 15    | 85    |

**Flow Rate:** 0.27 mL/min

**Col. Temp.:** ambient


**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 220 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Widespore C4 4 x 2.0mm, 10/Pk Part No.: AJ0-4329

- Holder Part No.: KJ0-4282

14935

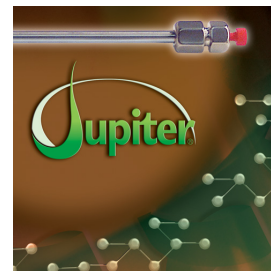

Products used in this application:

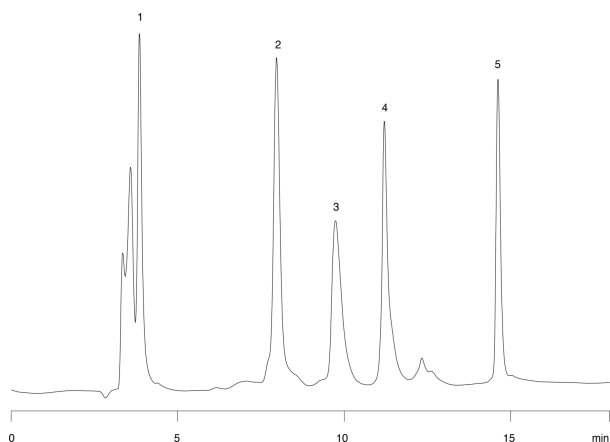

### ANALYTES:

- 1** Aprotinin  
Retention Time: 3.85 min
- 2** Ribonuclease  
Retention Time: 7.96 min
- 3** Fibrinogen  
Retention Time: 9.73 min
- 4** Acid glycoprotein  
Retention Time: 11.2 min
- 5** Leptin  
Retention Time: 14.61 min

