

## Triazine Herbicides on ZB-5MS

**Column:** Zebtron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea

**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane

**Dimensions:** 30 meters x 0.25 mm x 0.25 μm

**Order No:** 7HG-G010-11

**Oven Profile:** 140C to 220°C at 4°C/min

**Carrier Gas:** Constant Flow Helium, 1.2 mL/min


**Injection:** split 20:1 1 μL @ 250°C

**Detection:** Mass Selective (MSD) (240°C)

14954



Products used in this application:

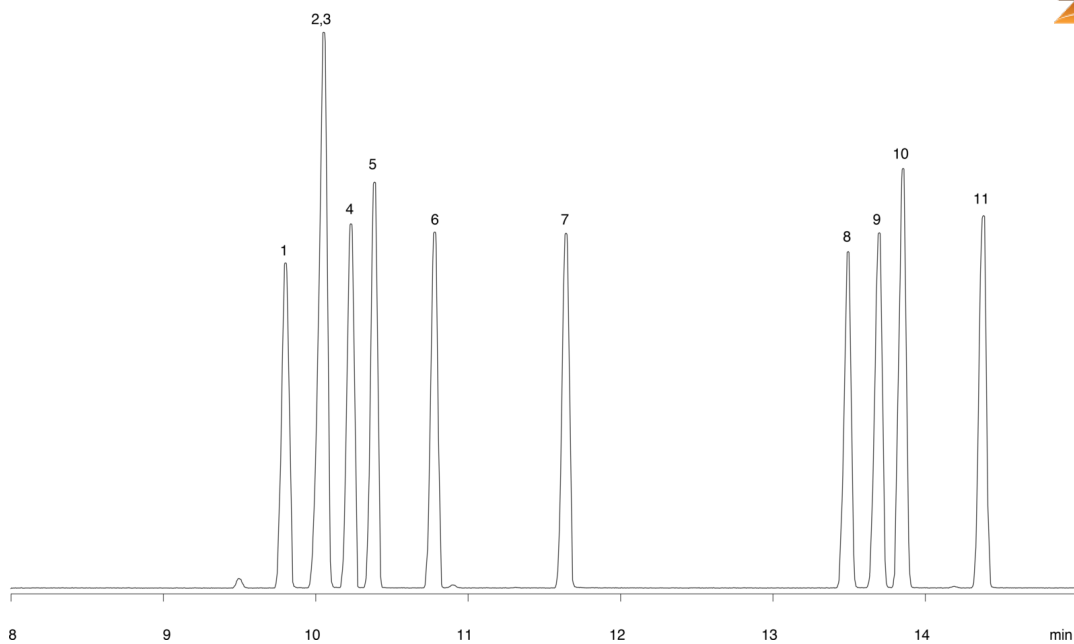

## Triazine Herbicides on ZB-5MS

### ANALYTES:

- 1 Atraton
- 2 Simazine
- 3 Prometon
- 4 Atrazine
- 5 Propazine
- 6 Terbutylazine
- 7 Secbumetone
- 8 Simetryn
- 9 Ametryn
- 10 Prometryn
- 11 Terbutryn

