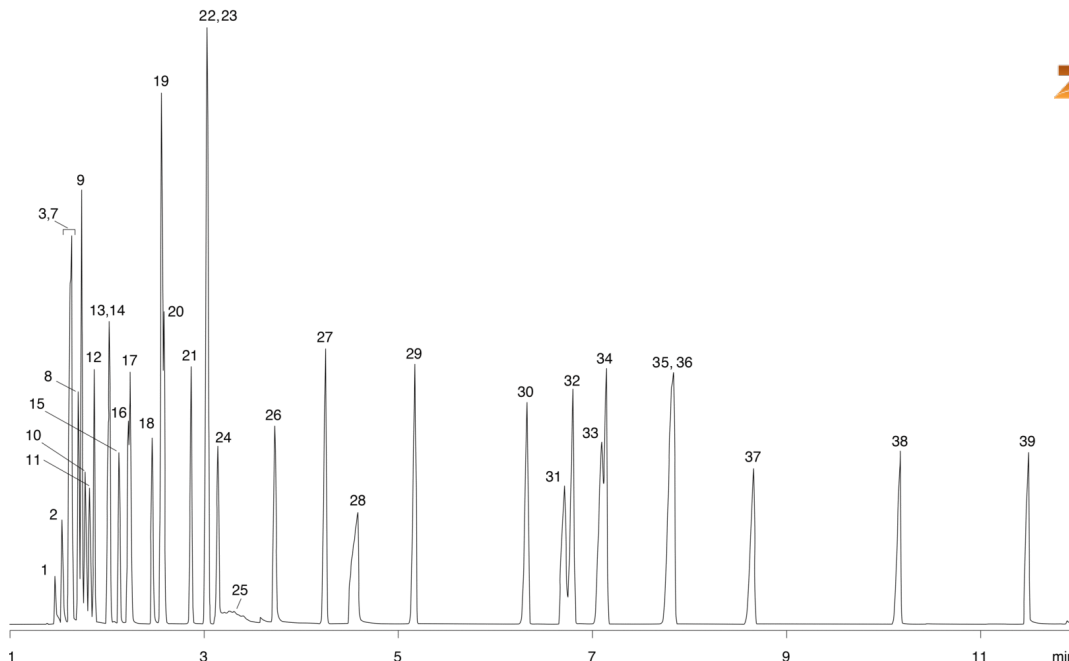


Solvents (mix 1) on ZB-5MS**Column:** Zebron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G010-11**Oven Profile:** 40°C for 2 min to 70°C at 4°C/min to 280°C at 25°C/min**Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** split 100:1 0.2 µL @ 250°C**Detection:** Mass Selective (MSD) (240°C)

15418

**Products used in this application:**

Solvents (mix 1) on ZB-5MS**ANALYTES:**

- 1** Methanol
- 2** Ethanol
- 3** Acetonitrile
- 4** Acetone
- 5** Pentane
- 6** Isopropanol
- 7** Ethyl ether
- 8** tert-Butanol
- 9** Methylene chloride
- 10** Allyl alcohol
- 11** n-Propanol
- 12** Methyl tert-butyl ether
- 13** 2-Butanone
- 14** Hexane
- 15** Ethyl acetate
- 16** Chloroform
- 17** Tetrahydrofuran
- 18** 1,2-Dimethoxy ethane
- 19** Benzene
- 20** Cyclohexane
- 21** isooctane
- 22** Heptane
- 23** Trichloroethylene
- 24** 1,4-Dioxane
- 25** Ethylene glycol
- 26** Pyridine
- 27** Toluene
- 28** N,N-Dimethylformamide
- 29** Tetrachloroethylene
- 30** Chlorobenzene
- 31** 3-Methylpyridine
- 32** Ethylbenzene
- 33** m-Xylene
- 34** p-Xylene
- 35** Styrene
- 36** o-Xylene
- 37** Anisole
- 38** Benzaldehyde
- 39** N-Methylpyrrolidinone

