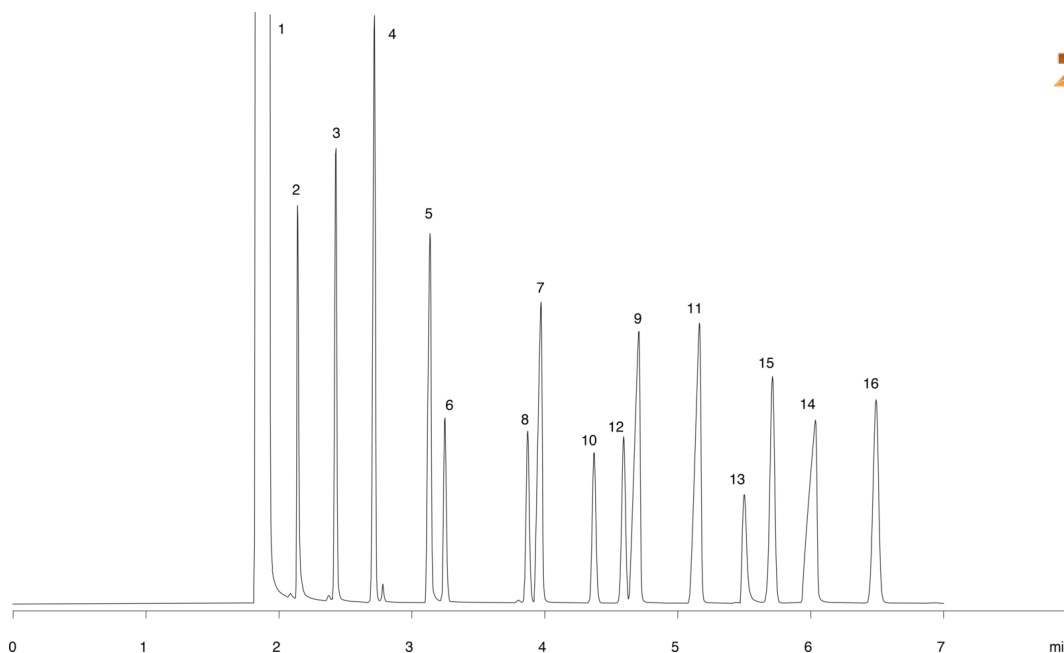


Petroleum Oxygenates on ZB-5MS**Column:** Zebtron[™] ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.50 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.5 µm**Order No:** 7HG-G010-17**Oven Profile:** 30°C Isothermal**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** split 100:1 0.2 µL @ 200°C**Detection:** Flame Ionization (FID) (260°C)

15434

**Products used in this application:**

Petroleum Oxygenates on ZB-5MS

ANALYTES:

- 1 Methanol
- 2 Ethanol
- 3 2-Propanol
- 4 tert-Butanol
- 5 n-Propanol
- 6 Methyl tert-butyl ether
- 7 Isopropyl ether
- 8 sec-Butanol
- 9 tert-Butyl ethyl ether
- 10 Methylcyclopentane
- 11 Isobutanol
- 12 tert-Pentanol
- 13 1,2-Dimethoxyethane
- 14 Benzene
- 15 n-Butanol
- 16 tert-Amyl methyl ether

