


Pharmaceutical Mix on Synergi 2.5µm Fusion-RP 20 x 2.0mm Mini-Column**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1% Formic acid in Water**Eluent B:** 0.1% Formic acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	92	8
	2	1	5	95
	3	1.01	92	8
	4	2	92	8

Flow Rate: 0.6 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**ANALYTES:**

- 1** Acetaminophen
Retention Time: 0.38 min
- 2** Caffeine
Retention Time: 0.634 min
- 3** Propanolol
Retention Time: 0.815 min
- 4** Imipramine
Retention Time: 0.884 min
- 5** Prednisolone
Retention Time: 1.037 min
- 6** Naproxen
Retention Time: 1.085 min
- 7** Indomethacine
Retention Time: 1.174 min

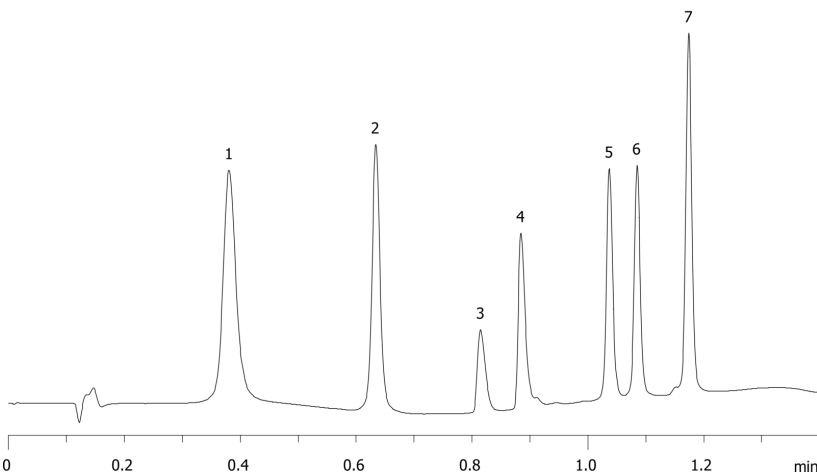