


Pharmaceutical Compounds on Synergi Fusion-RP 2.5µm 20 x 2.0mm MercuryMS Cartridge

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1% Formic Acid in Water
Eluent B: 0.1% Formic acid in Acetonitrile
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	92	8
2	1	5	95
3	1.01	92	8
4	2	92	8

Flow Rate: 0.6 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
 15516

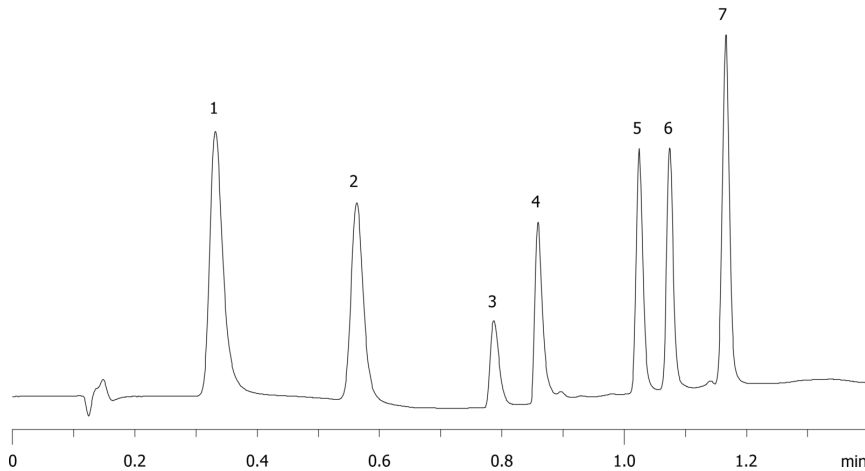

ANALYTES:

- 1** Acetaminophen
Retention Time: 0.332 min
- 2** Caffeine
Retention Time: 0.563 min
- 3** Propanolol
Retention Time: 0.787 min
- 4** Imipramine
Retention Time: 0.859 min
- 5** Prednisolone
Retention Time: 1.024 min
- 6** Naproxen
Retention Time: 1.074 min
- 7** Indomethacine
Retention Time: 1.165 min

