

## Catecholamines metabolites using Synergi Fusion-RP

**Column:** Synergi™ 4 µm Fusion-RP 80 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4424-E0

**Elution Type:** Gradient

**Eluent A:** 20mM Ammonium formate, pH 3.5

**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	10	50	50

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Fusion-RP 4 x 3.0mm, 10/Pk Part No.: AJ0-7557

- Holder Part No.: KJ0-4282



15739

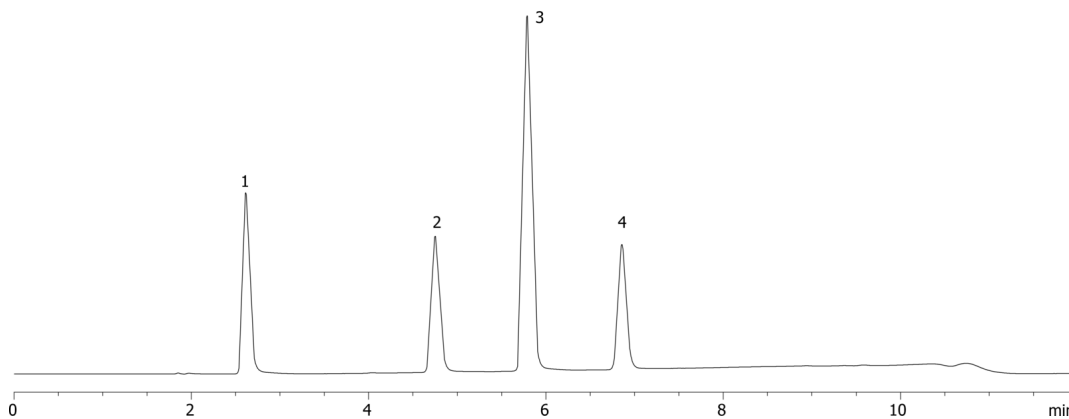

### ANALYTES:

- 1 Vanillylmandelic acid
- 2 3,4-Dihydroxyphenylacetic acid (DOPAC)
- 3 5-Hydroxyindoleacetic acid
- 4 Homovanillic acid

