

## Phenylacetic acids using Gemini 5u C6-Phenyl

**Column:** Gemini® 5 µm C6-Phenyl 110 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4444-E0

**Elution Type:** Isocratic

**Eluent A:** 20mM Potassium Phosphate pH 2.5

**Eluent B:** Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 78    | 22    |

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Gemini C6-Phenyl 4 x 3.0mm, 10/Pk Part No.: AJ0-7915


- Holder Part No.: KJ0-4282



Products used in this application:



15893

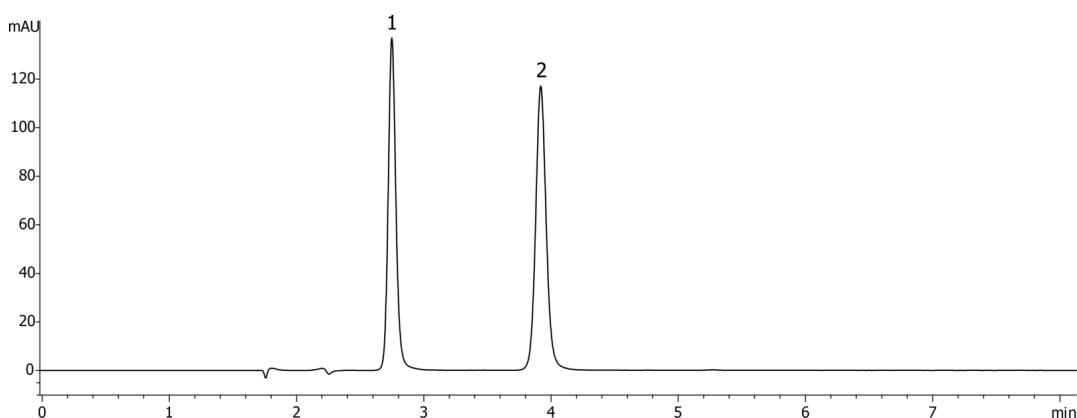

### ANALYTES:

- 1 3,4-Dihydroxy phenylacetic acid
- 2 4-Hydroxy-3-methoxyphenylacetic acid

