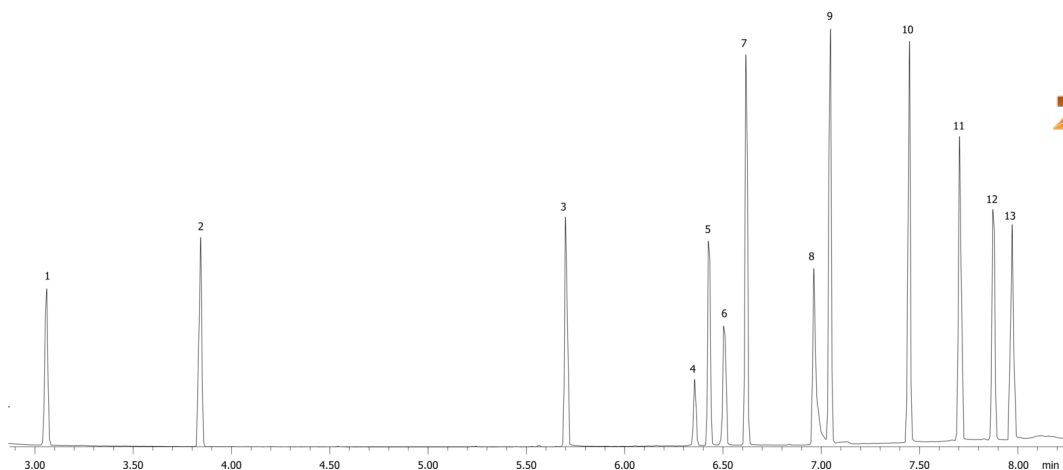


Grob Mix on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 45°C to 75°C at 6°C/min to 220°C at 30°C/min for 3 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 40:1 0.2 μ L @ 250°C**Detection:** Mass Selective (MSD) (250°C)

15909

**Products used in this application:**

Grob Mix on ZB-WAXplus

ANALYTES:

- 1 Decane
- 2 Undecane
- 3 Nonanal
- 4 2,3-Butanediol
- 5 1-Octanol
- 6 2,3-Butanediol
- 7 Methyl decanoate
- 8 Dicyclohexylamine
- 9 Methyl undecanoate
- 10 Methyl dodecanoate
- 11 2,6-Dimethylaniline
- 12 2,6-Dimethylphenol
- 13 2-Ethylhexanoic acid

