


Ethylene and Propylene Glycols on ZB-WAXplus**Column:** Zebron[™] ZB-WAXplus[™], GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G013-11**Oven Profile:** 80°C to 200°C at 8°C/min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Splitless :1 0.2 μ L @ 140°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Analytes are 100ppm in HPLC water**Products used in this application:**

15911

**ANALYTES:**

- 1 Propylene glycol
- 2 Ethylene glycol
- 3 HPLC water contaminant

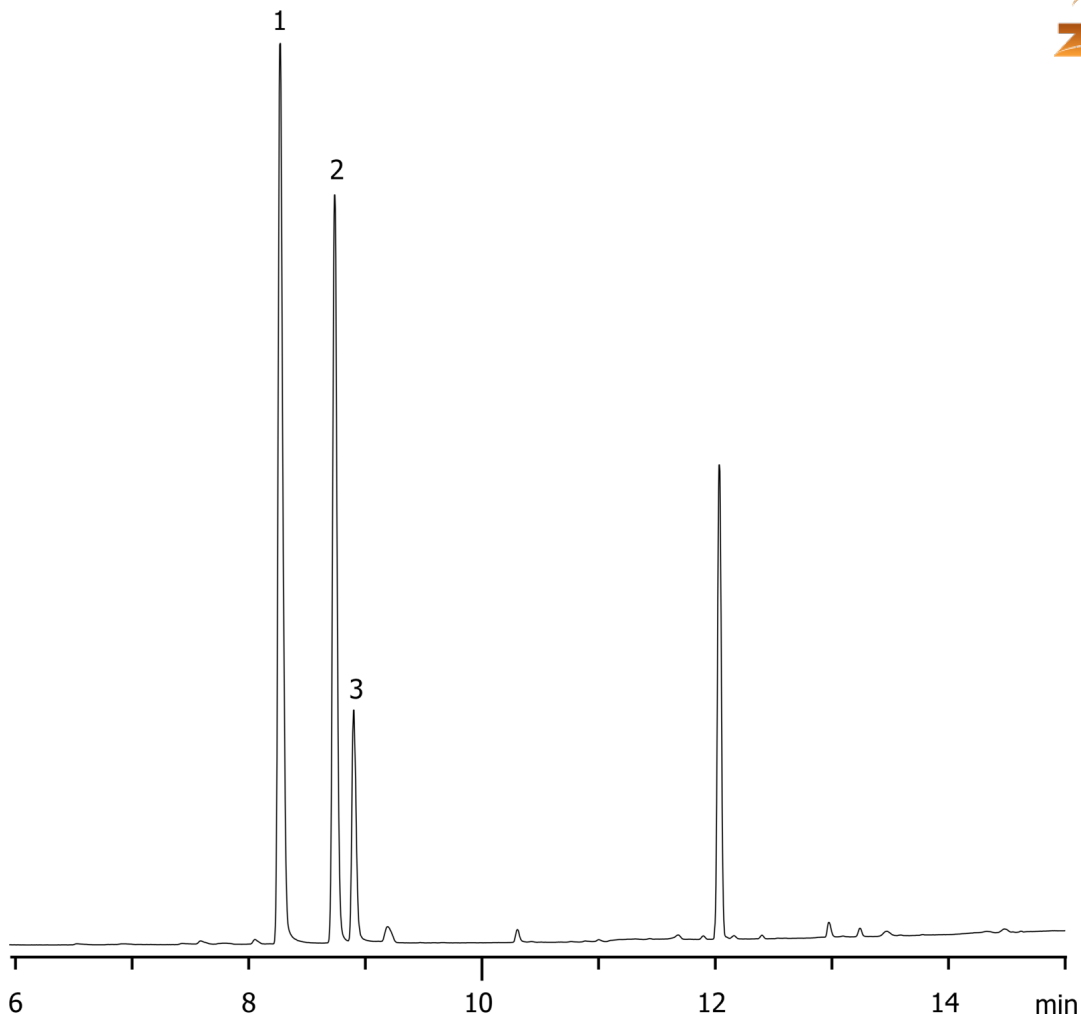