


Solvents on ZB-WAXplus 30x0.25x0.25**Column:** Zebtron™ ZB-WAXplus™, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G013-11**Oven Profile:** 30°C for 1min to 70°C at 14°C/min to 220°C at 25°C/min for 3 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 50:1 0.2 µL @ 220°C**Detection:** Mass Selective (MSD) (250°C)**Analyst Note:** All analytes are 85 ppm in DMA**Products used in this application:**

Solvents on ZB-WAXplus 30x0.25x0.25**ANALYTES:**

1 air	41 Butanol
2 Pentane	42 Nitromethane
3 Hexane	43 Cumene
4 Ethyl ether	44 2-Methoxy ethanol
5 2-Methoxy-2-methyl-2-propane	45 Pyridine
6 Heptane	46 o-Xylene
7 Cyclohexane	47 3-Methyl butanol
8 1,1-Dichloroethene	48 Chlorobenzene
9 Methylcyclohexane	49 n-Pentanol
10 Acetone	50 DMF
11 Formic acid ethyl ester	51 Anisole
12 Acetic acid methyl ester (methyl acetate)	52 Tetralin
13 Tetrahydrofuran	53 Dimethyl sulfoxide
14 Carbon tetrachloride	54 Ethylene glycol
15 1,1,1-Trichloroethane	55 Methylacetamide
16 Ethyl acetate	56 N-Methylpyrrolidinone
17 Acetic acid-1-methyl ethyl ester	57 Formamide
18 2-Butanone	
19 1,2-Dimethoxyethane	
20 Dichloromethane	
21 Isopropanol	
22 Benzene	
23 Propyl acetate	
24 1,2-Dichloroethene	
25 Trichloroethylene	
26 Acetonitrile	
27 Methylisobutylketone (MIBK)	
28 Isobutyl acetate	
29 Chloroform	
30 Toluene	
31 n-Propanol	
32 Water	
33 1,4-Dioxane	
34 1,2-Dichloroethane	
35 Butyl acetate	
36 2-Hexanone	
37 2-Methylpropanol	
38 Ethylbenzene	
39 p-Xylene	
40 m-Xylene	

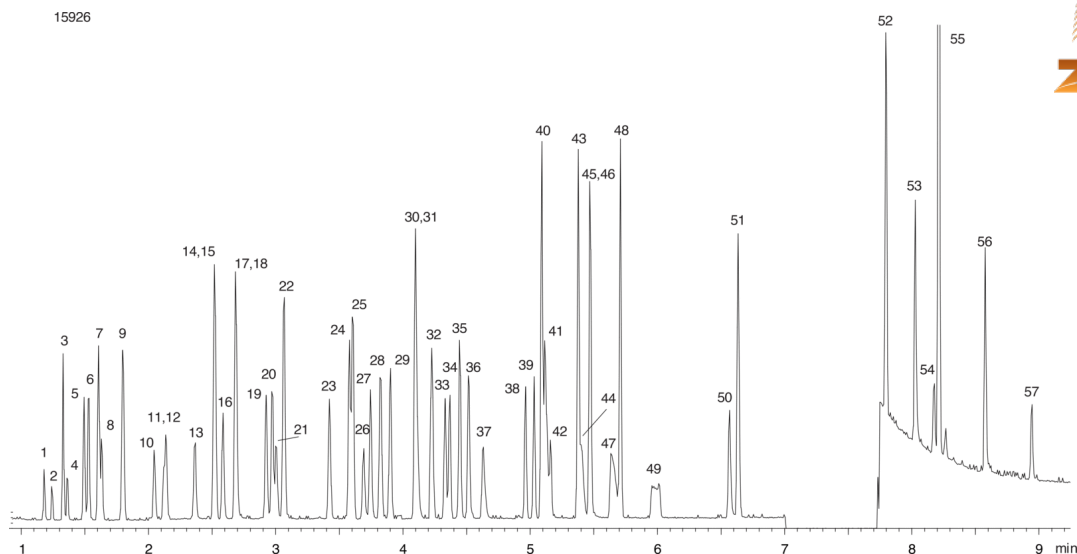