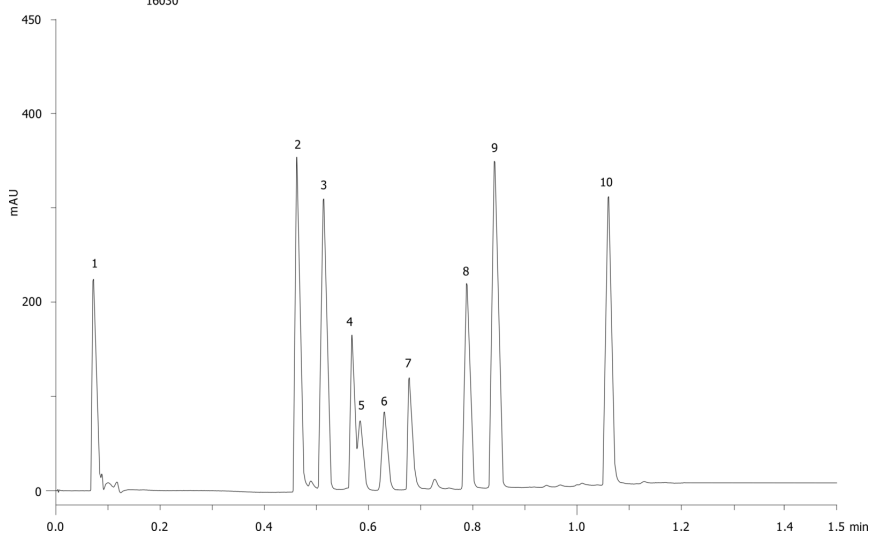


MercuryMS - Synergi 2.5µm Fusion-RP 20 x 4.0mm for Acid, Base, Neutral

Column: N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1	5	95
	3	1.01	95	5
	4	1.2	95	5

Flow Rate: 2.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

MercuryMS - Synergi 2.5 μ m Fusion-RP 20 x 4.0mm for Acid, Base, Neutral

ANALYTES:

- 1 Pyridine
- 2 Quinidine
- 3 Sulfathiazole
- 4 Triprolidine
- 5 Benzyl alcohol
- 6 Phenol
- 7 Nortriptyline
- 8 3-Methyl-4-nitrobenzoic acid
- 9 5-Methyl-salicylaldehyde
- 10 Hexaphenone

