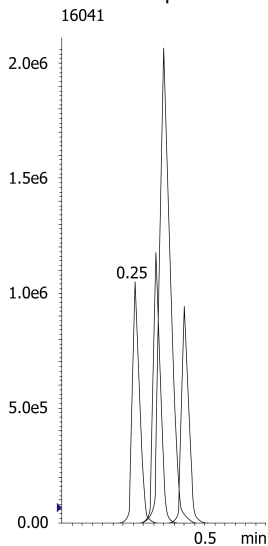


MercuryMS - Synergi 2.5µm Fusion-RP 20 x 2.0mm for Tricyclics in Plasma by LC/MS/MS**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic Acid in Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	28	72

Flow Rate: 0.6 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**ANALYTES:****1** Doxepin

Retention Time: 0.42 min

2 Desipramine

Retention Time: 0.56 min

3 Imipramine

Retention Time: 0.63 min

4 Nortriptyline

Retention Time: 0.67 min

5 Trimipramine

Retention Time: 0.85 min

6 Amitriptyline

Retention Time: 0.76 min

