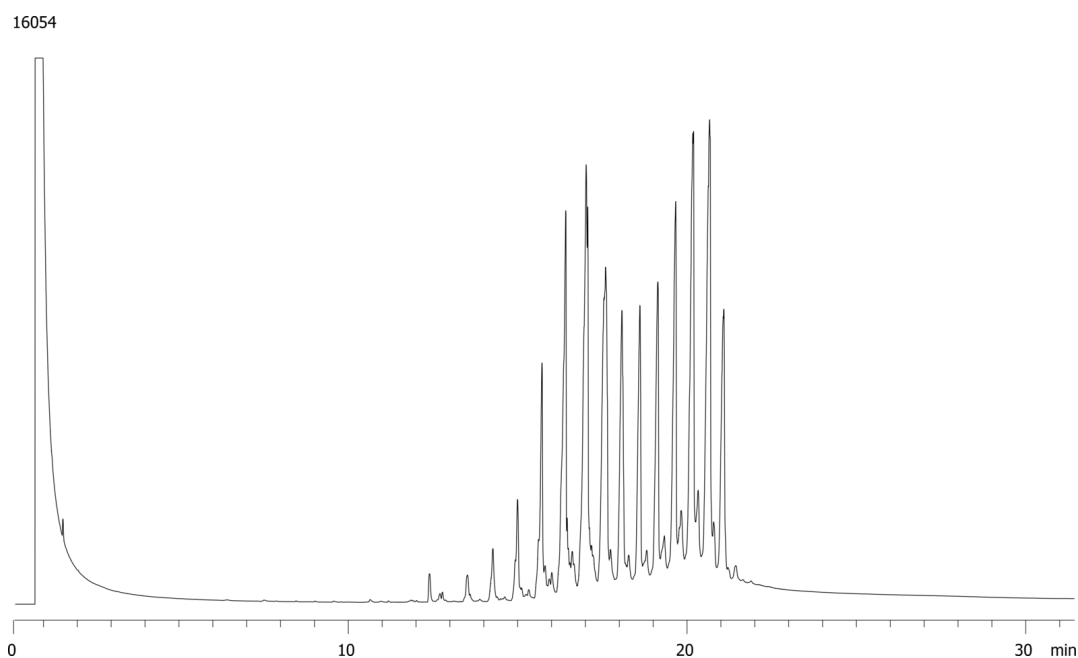


Butter Triglycerides on ZB-5HT 15m x 0.32mm x 0.10µm**Column:** Zebtron[™] ZB-5HT, GC Cap. Column 15 m x 0.32 mm x 0.10 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.1 µm**Order No:** 7EM-G015-02**Oven Profile:** 100°C to 400°C at 14°C/min for 10 min**Carrier Gas:** Constant Flow Helium, 1.8 mL/min**Injection:** On-Column :1 2 µL @ 103°C**Detection:** Flame Ionization (FID) (400°C)**Analyst Note:** Recommended Accessories:
Recommended Liner: Zebtron PLUS Liner for Agilent & Thermo
4mm ID Direct Connect

Liner Part No.: AG2-0A50-05

**Products used in this application:****ANALYTES:****1** butter