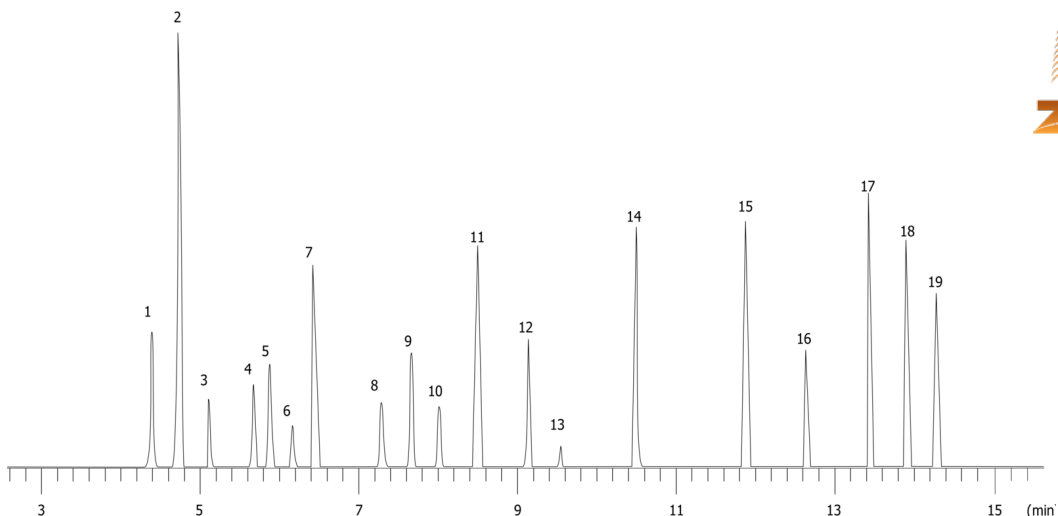


Solvents on ZB-FFAP**Column:** ZebronTM ZB-FFAP, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Nitroterephthalic Acid Modified Polyethylene Glycol**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G009-11**Oven Profile:** 60 °C for 8 min to 150 °C at 15 °C/min for 4 min**Carrier Gas:** Constant Pressure Helium, 3.4 mL/min**Injection:** Split 11.8:1 1 μ L @ 225°C**Detection:** Flame Ionization (FID) (250°C)

16215

**Products used in this application:**

Solvents on ZB-FFAP

ANALYTES:

- 1 n-Hexane
- 2 Carbon disulfide
- 3 Acetone
- 4 Ethyl acetate
- 5 Methyl ethyl ketone
- 6 Dichloromethane
- 7 Benzene
- 8 Trichloroethylene
- 9 Methylisobutylketone
- 10 Perchloroethylene
- 11 Toluene
- 12 n-Butyl acetate
- 13 Undecane
- 14 Ethylbenzene
- 15 o-Xylene
- 16 PGMEAC
- 17 Styrene
- 18 1,2,4-Trimethylbenzene
- 19 Cyclohexane

