


TMS Sugars on ZB-MR1**Column:** Zebtron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 200 °C for 6 min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Split 100:1 1 μ L @ 225°C**Detection:** Flame Ionization (FID) (250°C)**Analyst Note:** Sugar mix concentration 176 ppm. Sample was blown dry and then reconstituted in 100 μ L of 1:1**Products used in this application:**

16248

**ANALYTES:**

- 1 D-(-)-Arabinose
- 2 d-Ribose
- 3 d-Xylose
- 4 D-Mannose
- 5 D-(+)-Galactose
- 6 alpha-d-Glucose
- 7 beta-D-(+)-Glucose

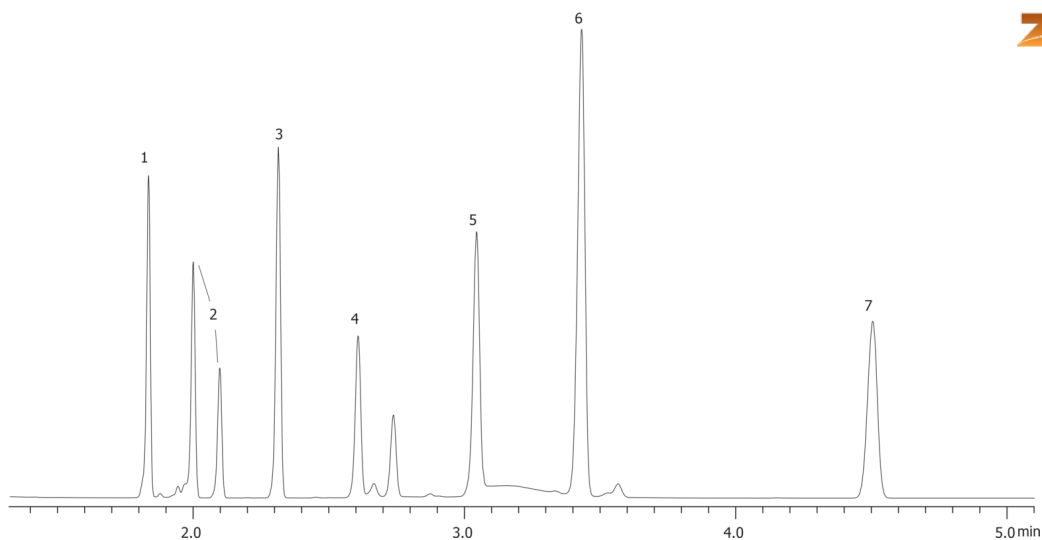