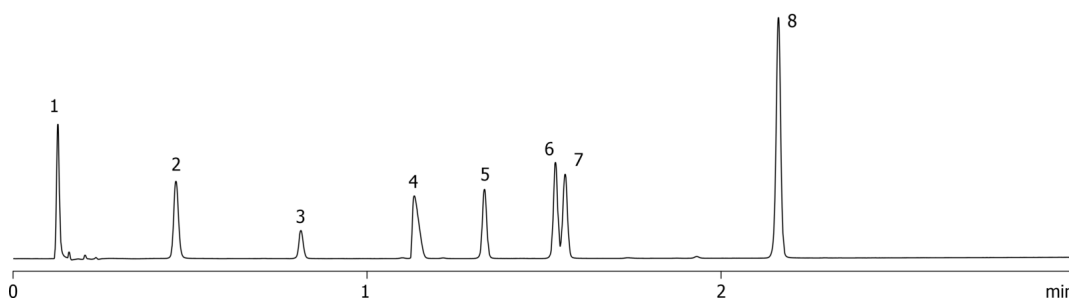


Acid, Base, Neutral Test Mix on Synergi 2.5 μ m Fusion-RP-HST**Column:** SynergiTM 2.5 μ m Fusion-RP 100 Å, LC Column 50 x 2 mm, Ea**Dimensions:** 50 x 2 mm ID**Order No:** 00B-4423-B0**Elution Type:** Gradient**Eluent A:** 0.1 % Formic acid in Water**Eluent B:** 0.1 % Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2.9	5	95

Flow Rate: 1.1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)**Detector Info:** Agilent 1200 SL

16265

**ANALYTES:**

- 1 Pyridine
- 2 Acetaminophen
- 3 Benzyl alcohol
- 4 Nortriptyline
- 5 3-Methyl-4-nitrobenzoic acid
- 6 4-Chlorocinnamic acid
- 7 3-Hydroxy-3-methylbenzaldehyde
- 8 Hexanophenone

