

## Comparison of Water-Soluble Vitamins: Luna 5µm HILIC 90/10 ACN/Buffer, pH 3.2

**Column:** Luna® 5 µm HILIC 200 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4450-E0

**Elution Type:** Isocratic

**Eluent A:** Acetonitrile

**Eluent B:** 100mM Ammonium Formate pH 3.2

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10

**Flow Rate:** 1 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 260 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, HILIC 4 x 3.0 mm, 10/Pk Part No.: AJ0-8329

- Holder Part No.: KJ0-4282



Products used in this application:

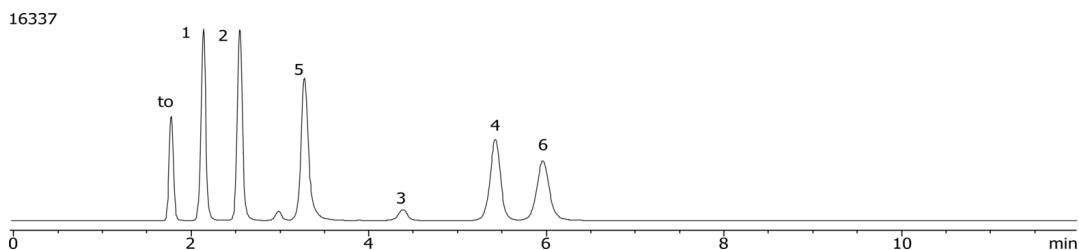

### ANALYTES:

- 1 PABA
- 2 Nicotinamide
- 3 Riboflavin
- 4 Nicotinic acid
- 5 Pyridoxine
- 6 Thiamine

