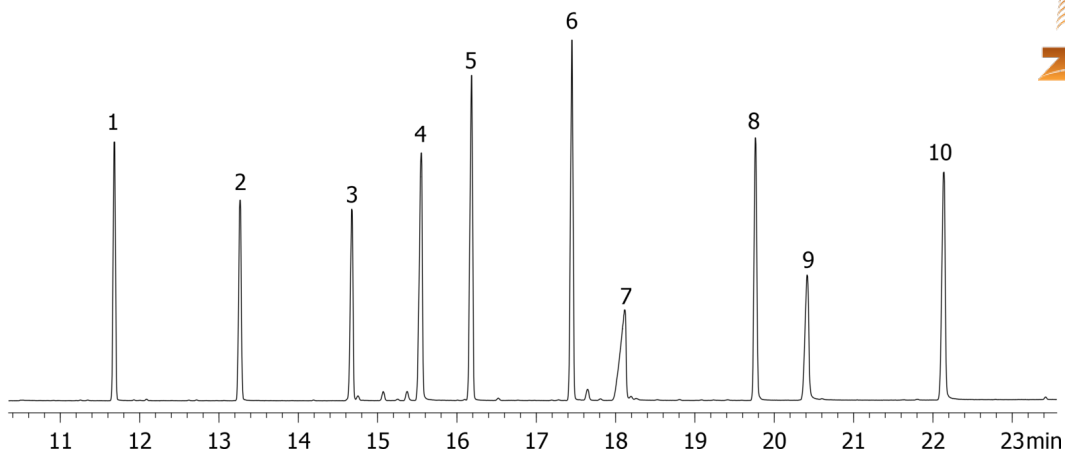


Nitrogen and Phosphorous Pesticides on ZB-5MS**Column:** ZebtronTM ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G010-11**Oven Profile:** 60 °C for 0.5 min to 320°C @ 10 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 0.9 mL/min**Injection:** Splitless :1 1 μ L @ 260°C**Detection:** Mass Selective (MSD) (180°C)

16368

**Products used in this application:**

Nitrogen and Phosphorous Pesticides on ZB-5MS

ANALYTES:

- 1 Butylate
- 2 Molinate
- 3 Chlorpropham
- 4 Atraton
- 5 Pronamide
- 6 Alachlor
- 7 Bromacil
- 8 Tetrachlorvinphos
- 9 Tricyclazole
- 10 Hexazinone

