

## Nitrogen and Phosphorous Pesticides on ZB-5MS

**Column:** Zebtron™ ZB-5MS, GC Cap. Column 30 m x 0.25 mm x 0.25 μm, Ea

**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane

**Dimensions:** 30 meters x 0.25 mm x 0.25 μm

**Order No:** 7HG-G010-11

**Oven Profile:** 60 °C for 0.5 min to 320°C @ 10 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 0.9 mL/min


**Injection:** Splitless :1 1 μL @ 260°C

**Detection:** Mass Selective (MSD) (180°C)

16368



Products used in this application:

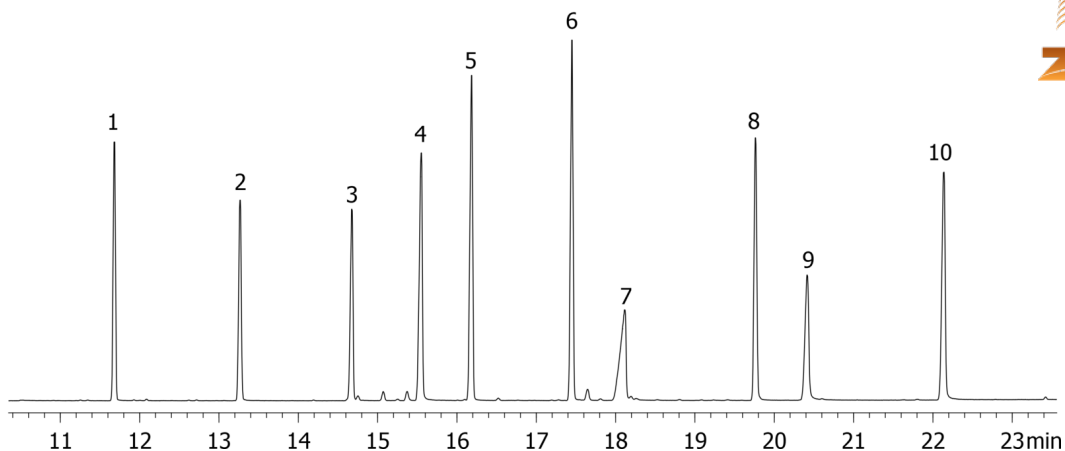

## Nitrogen and Phosphorous Pesticides on ZB-5MS

### ANALYTES:

- 1 Butylate
- 2 Molinate
- 3 Chlorpropham
- 4 Atraton
- 5 Pronamide
- 6 Alachlor
- 7 Bromacil
- 8 Tetrachlorvinphos
- 9 Tricyclazole
- 10 Hexazinone

