

Nitrogen and Phosphorous Pesticides on ZB-MR2

Column: Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μm, Ea

Phase: Unknown or Proprietary

Dimensions: 30 meters x 0.25 mm x 0.2 μm

Order No: 7HG-G017-10

Oven Profile: 60 °C for 0.5 min to 120 °C @ 10 °C/min to 240 °C @ 4 °C/min to 320 °C @ 15 °C/min for 1 min

Carrier Gas: Constant Flow Helium, 0.9 mL/min


Injection: Splitless :1 1 μL @ 260°C

Detection: Mass Selective (MSD) (180°C)

16383



Products used in this application:

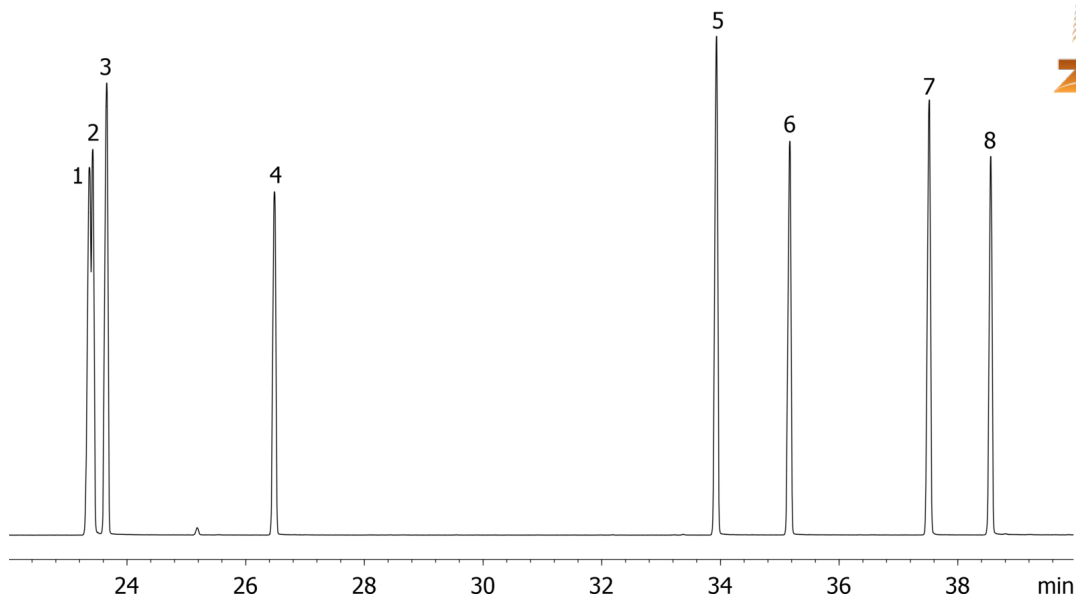

ANALYTES:

- 1 Propachlor
- 2 Trifluralin
- 3 Benefin (Benfluraline)
- 4 Profluralin
- 5 Isopropalin
- 6 Pendimethalin
- 7 Oxadiazon
- 8 Oxyfluorfen

