


Rose Oil on ZB-XLB**Column:** Zebron[™] ZB-XLB, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary XLB Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G019-11**Oven Profile:** 50 °C to 340 °C at 10 °C/min for 10 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Split 25:1 1 μ L @ 250°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** 1% Rose oil in ethyl acetate**Products used in this application:**

Rose Oil on ZB-XLB

ANALYTES:

- 1 α -Pinene
- 2 2-Ethylhexanol
- 3 Benzyl alcohol
- 4 Terpineol
- 5 Linalool
- 6 Phenylethyl alcohol
- 7 β -Citronellol
- 8 cis-Geraniol
- 9 trans-Geraniol
- 10 β -Citral
- 11 Geranic acid
- 12 Eugenol
- 13 Geranyl acetate
- 14 Eugenol methyl ether
- 15 Caryophyllene
- 16 α -Guaiene
- 17 α -Caryophyllene
- 18 Germacrene D
- 19 α -Bulnesene
- 20 Heptadecane (C17)
- 21 Farnesol
- 22 1-Nonadecene
- 23 Nonadecane
- 24 Vitamin E
- 25 α -Amirin

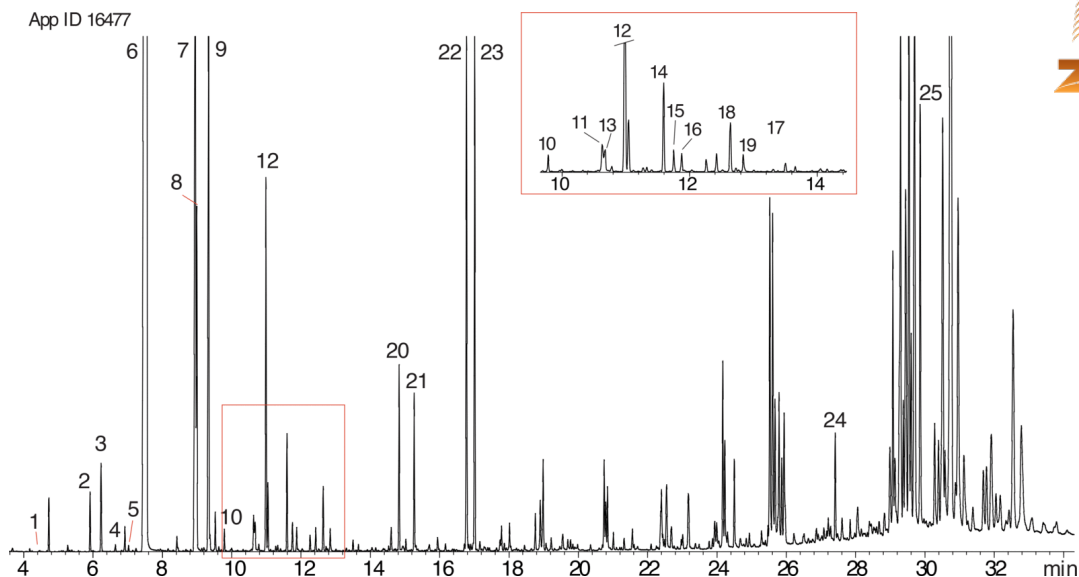