

Polar bases (beta blockers) at high pH (ammonium bicarbonate pH 10.5) on Gemini-NX 5μ

Column: Gemini® 5 μm NX-C18 110 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4454-E0

Elution Type: Gradient

Eluent A: 10mM Ammonium Bicarbonate pH 10.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	15	70	30
	3	20	50	50
	4	25	50	50

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

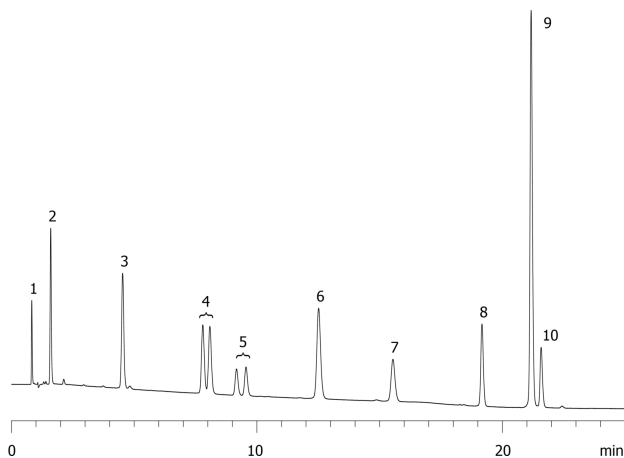

Polar bases (beta blockers) at high pH (ammonium bicarbonate pH 10.5) on Gemini-NX 5 μ

ANALYTES:

- 1** Bisoprolol
- 2** Sotalol
- 3** Atenolol
- 4** Labetalol
- 5** Labetalol
- 6** Nadolol
- 7** Nadolol
- 8** Pindolol
- 9** Metoprolol
- 10** Bisoprolol
- 11** Propranolol
- 12** Alprenolol

