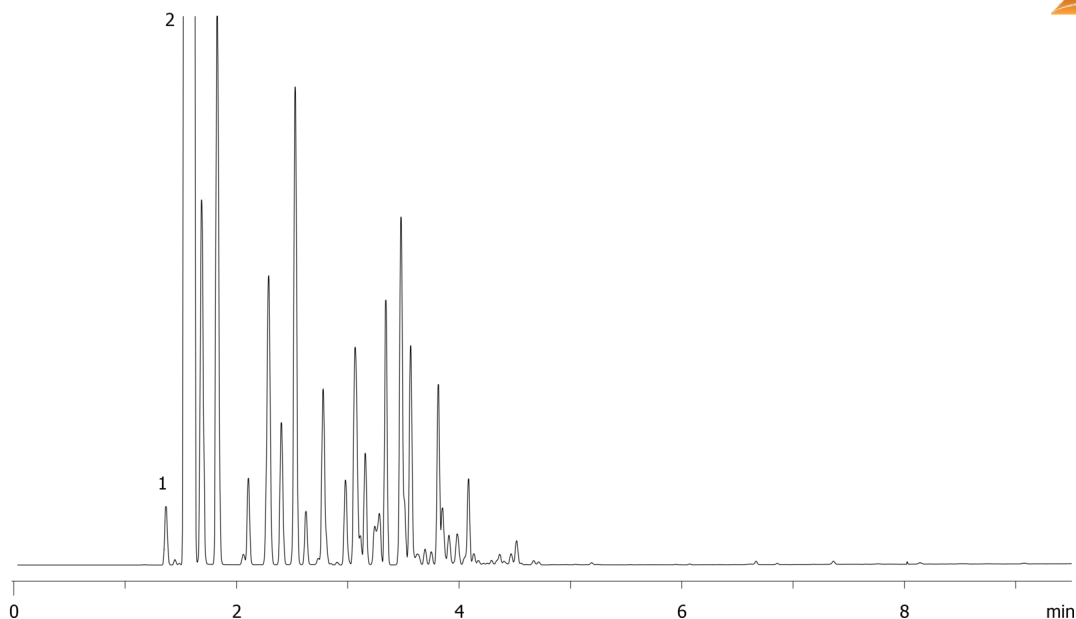


**ASTM D5501: Determination of Denatured Bioethanol on ZB-Bioethanol****Column:** Zebron<sup>™</sup> ZB-Bioethanol, GC Cap. Column 15 m x 0.25 mm x 1.00  $\mu$ m, Ea**Phase:** Unknown or Proprietary**Dimensions:** 15 meters x 0.25 mm x 1  $\mu$ m**Order No:** 7EG-G020-22**Oven Profile:** 55 °C for 1.7 min to 260 °C @ 40 °C/min for 2.67 min.**Carrier Gas:** Constant Flow Hydrogen, 25 cm/sec**Injection:** Split 50:1 1  $\mu$ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Instrumentation: Shimadzu GC-2010**Products used in this application:**

16668

**ANALYTES:****1** Methanol**2** Ethanol