

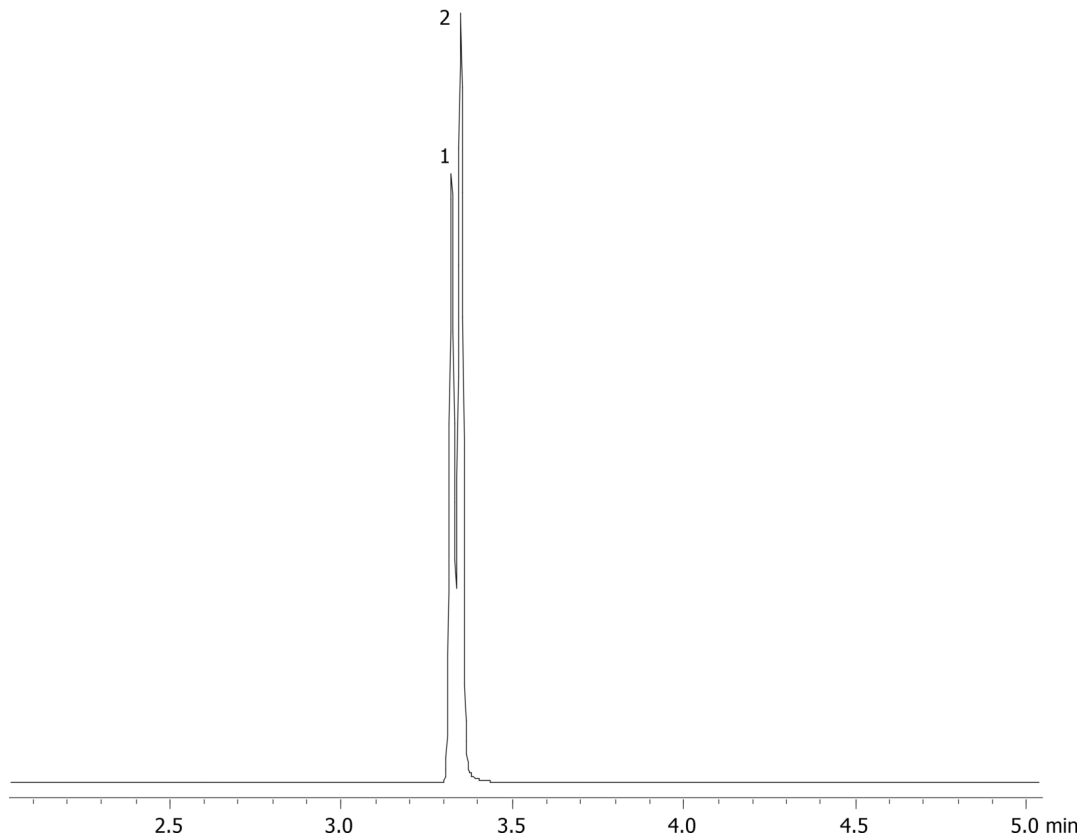
Methadone by MSD on a ZB-1MS

Column: Zebtron™ ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 µm, Ea
Phase: 100% Dimethylpolysiloxane
Dimensions: 10 meters x 0.18 mm x 0.18 µm
Order No: 7CD-G011-08
Oven Profile: 150 °C to 320 °C @ 20 °C/min

Carrier Gas: Constant Flow Helium, 1.2 mL/min
Injection: Split 10:1 1 µL @ 250°C
Detection: Mass Selective (MSD) (320°C)
Analyst Note: Analytes are 24ppm
Sim ions:
(Quant ion*) D9-Methadone (*78,303): Methadone (*72,294,223)

**Products used in this application:**

16692

**ANALYTES:**

- 1 Methadone-d9
- 2 Methadone

