


**Carboxy-THC by MSD on a ZB-5MS****Column:** Zebron<sup>™</sup> ZB-5MS, GC Cap. Column 10 m x 0.18 mm x 0.18  $\mu$ m, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18  $\mu$ m**Order No:** 7CD-G010-08**Oven Profile:** 200 °C to 300 °C @ 15 °C/ min**Carrier Gas:** Constant Flow He, 1.2 mL/min**Injection:** Split 10:1 1  $\mu$ L @ 250°C**Detection:** Mass Selective (MSD) (300°C)**Analyst Note:** Analytes are 1.5ppm derivatized with BSTFA  
Sim ions:  
(Quant ion\*) D3-Carboxy-thc (\*374,476,477): Carboxy-thc (\*371,473,474)**Products used in this application:****ANALYTES:**

- 1 Carboxy-THC-d3
- 2 Carboxy-THC

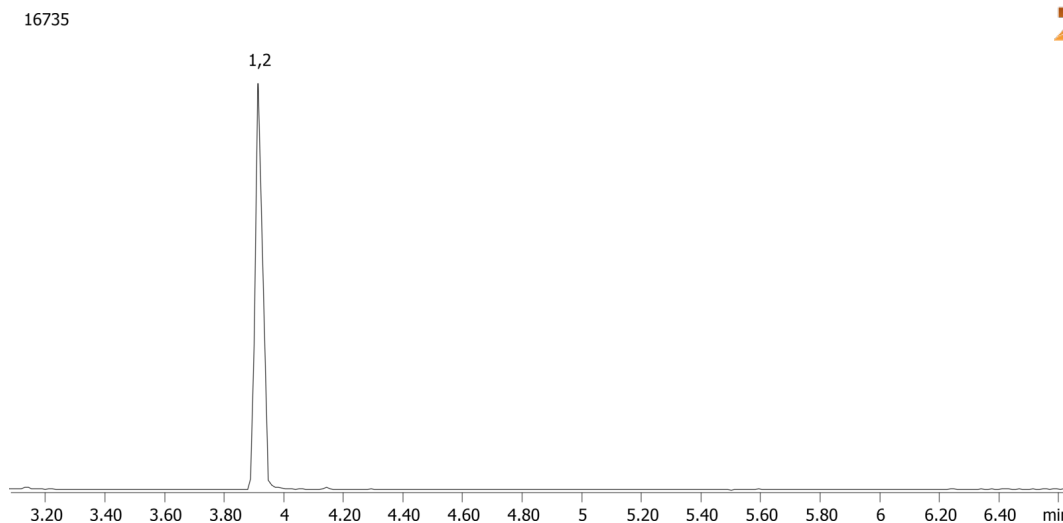