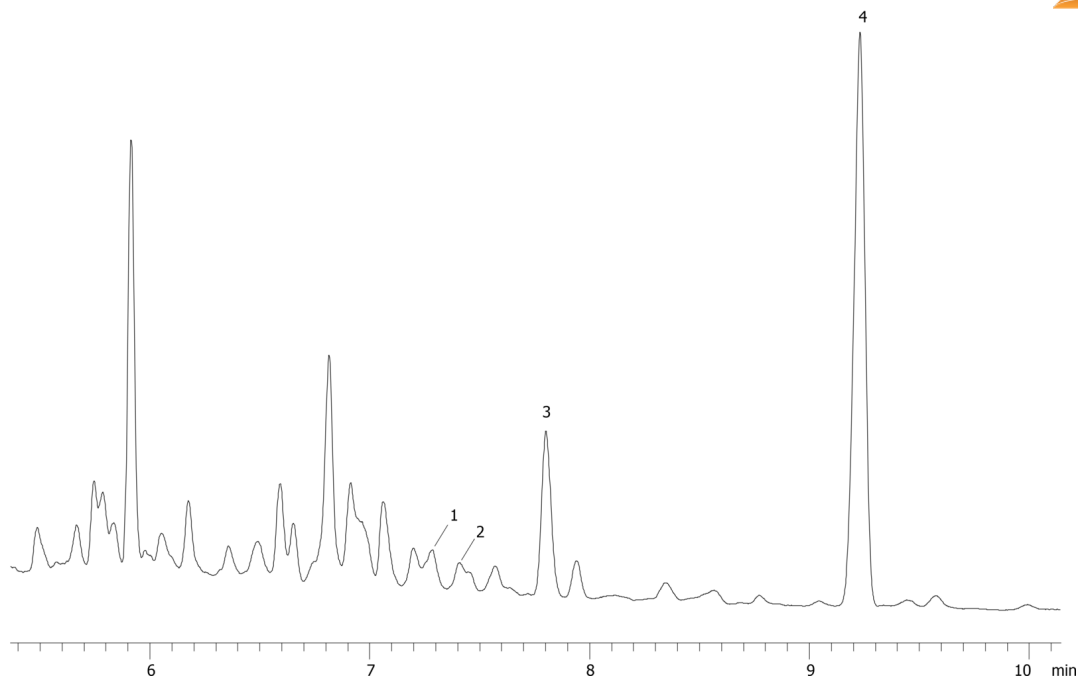


**Sterols from Peanut Butter by FID on ZB-MR1**

**Column:** Zebtron<sup>™</sup> ZB-MultiResidue<sup>™</sup>-1, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea  
**Phase:** Proprietary Pesticides Phase  
**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m  
**Order No:** 7HG-G016-11  
**Oven Profile:** 230 °C to 320 °C @ 15 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 2.5 mL/min  
**Injection:** Splitless :1 0.5  $\mu$ L @ 300°C  
**Detection:** Flame Ionization (FID) (300°C)  
**Analyst Note:** Sample prepared- Saponification, Solid Phase Extraction (SPE), and Derivatized.

16746

**Products used in this application:****ANALYTES:**

- 1 Campesterol
- 2 Stigmasterol
- 3 b-Sitosterol
- 4 Betulin

