


Sterols from Olive Oil by FID on ZB-MR1**Column:** Zebron[™] ZB-MultiResidue[™]-1, GC Cap. Column 30 m x 0.25 mm x 0.25 μ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 μ m**Order No:** 7HG-G016-11**Oven Profile:** 230 °C to 320 °C @ 25 °C/min for 5 min**Carrier Gas:** Constant Flow Helium, 2.5 mL/min**Injection:** Splitless :1 0.5 μ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sample Preparation- Saponification, pH Adjustment, Solid Phase Extraction (SPE), Reconstituted, and Derivatized.**Products used in this application:**

16747

**ANALYTES:**

- 1** Campesterol
- 2** Stigmasterol
- 3** b-Sitosterol
- 4** Betulin

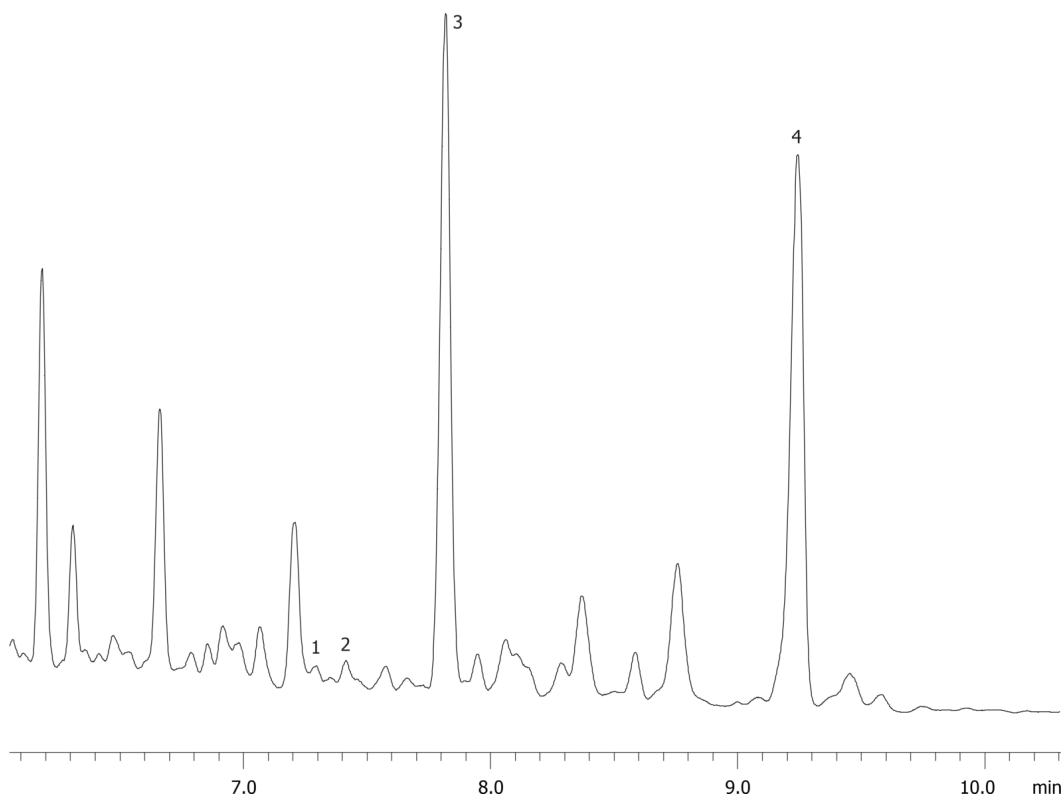