


Sterols from Olive Oil by FID on ZB-5HT**Column:** ZebtronTM ZB-5HT, GC Cap. Column 30 m x 0.25 mm x 0.10 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.1 μ m**Order No:** 7HG-G015-02**Oven Profile:** 220 °C to 350 °C @ 15 °C/min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5 μ L @ 350°C**Detection:** Flame Ionization (FID) (350°C)**Analyst Note:** Sample Prep: Saponification, Solid Phase Extraction (SPE), and Derivatization.
Recommended Accessories:
Recommended Liner: Zebtron PLUS Single Taper Z-Liner, 4 mm ID
Liner Part No.: AG2-0A13-05 (for Agilent systems)**Products used in this application:**

16753

**ANALYTES:**

- 1 Campesterol
- 2 Stigmasterol
- 3 b-Sitosterol
- 4 Betulin

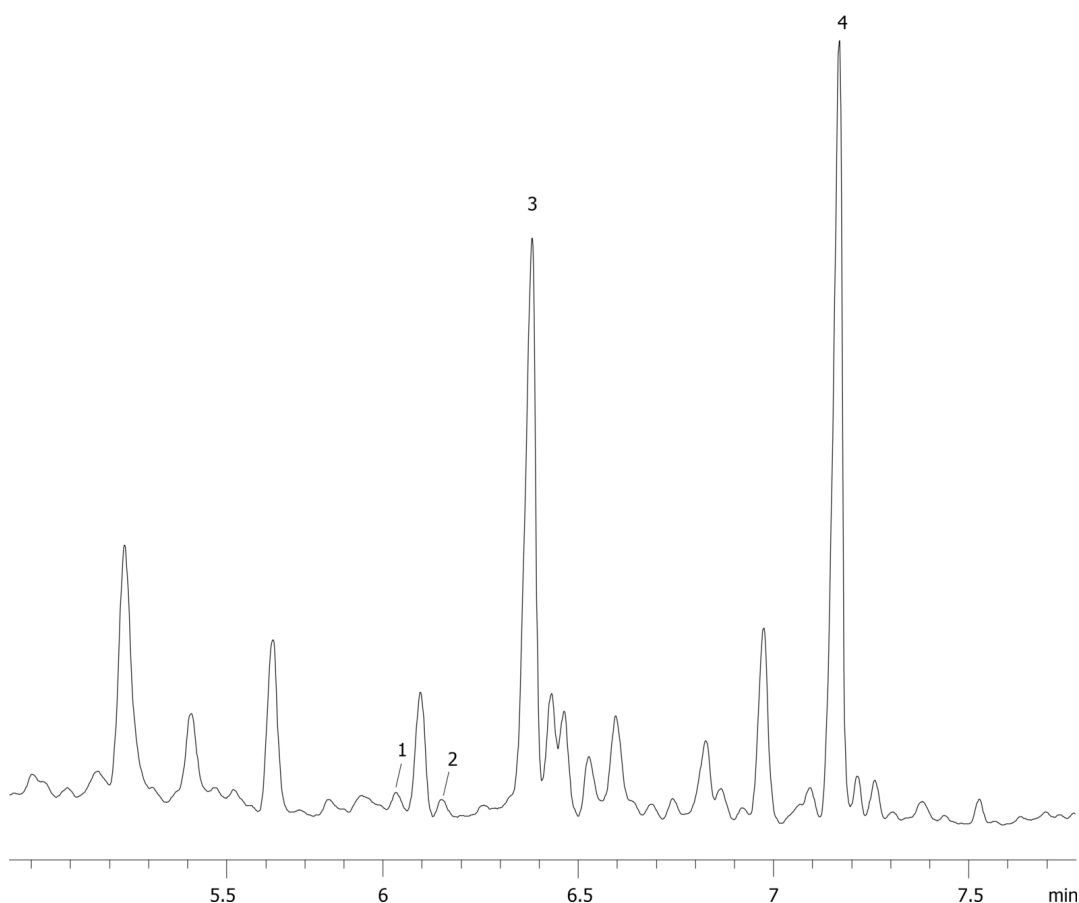