


**Sterols by FID on ZB-MR 2****Column:** Zebron<sup>TM</sup> ZB-MultiResidue<sup>TM</sup>-2, GC Cap. Column 30 m x 0.25 mm x 0.20  $\mu$ m, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.2  $\mu$ m**Order No:** 7HG-G017-10**Oven Profile:** 220 °C to 320 °C @ 13 °C/min for 5 min.**Carrier Gas:** Constant Flow Helium, 2 mL/min**Injection:** Splitless :1 0.5  $\mu$ L @ 300°C**Detection:** Flame Ionization (FID) (300°C)**Analyst Note:** Sterols derivatized via BSTFA:TMCS**Products used in this application:**

16754

**ANALYTES:**

- 1 Cholesterol
- 2 Campesterol
- 3 Stigmasterol
- 4 b-Sitosterol
- 5 Betulin

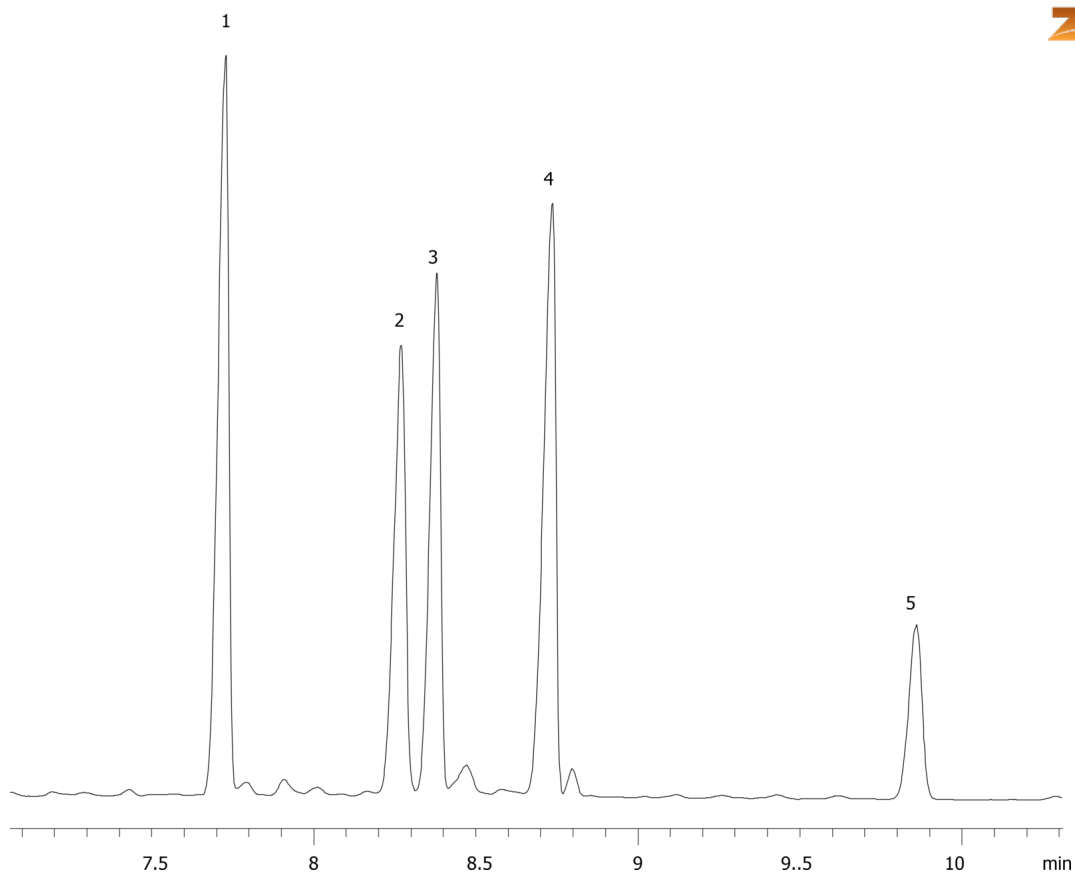