

Hydroxyzine on Lux Cellulose-1**Column:** Lux® 5 µm Cellulose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4459-E0**Elution Type:** Isocratic**Eluent A:** 0.1 % DEA in CH₃CN**Eluent B:** 0.1% DEA in IPA

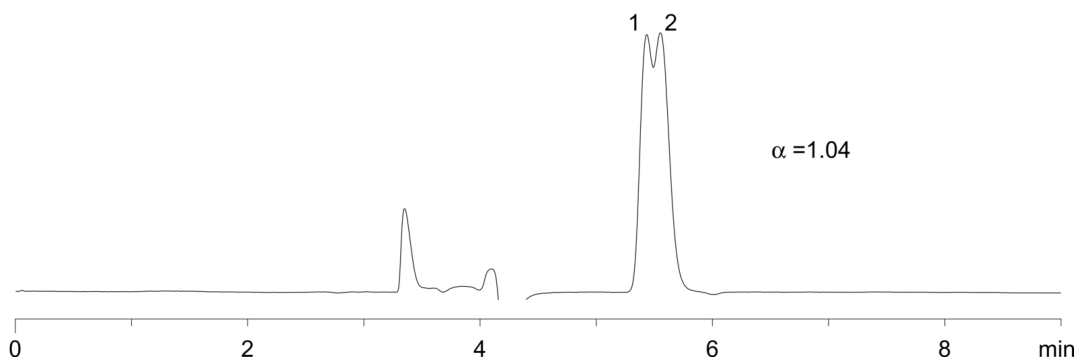
Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	10	100	0

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Lux Cellulose-1, 4 x 3.0mm, 10/Pk Part No.: AJ0-8403

- Holder Part No.: KJ0-4282

17495

**Products used in this application:****ANALYTES:****1** Hydroxyzine

Retention Time: 5.44 min

2 Hydroxyzine

Retention Time: 5.55 min

