

## Nitrogen and Phosphorous Pesticides (EPA 507) Mix # 2 on ZB-MR2

**Column:** Zebtron™ ZB-MultiResidue™-2, GC Cap. Column 30 m x 0.25 mm x 0.20 μm, Ea

**Phase:** Proprietary Pesticides Phase

**Dimensions:** 30 meters x 0.25 mm x 0.2 μm

**Order No:** 7HG-G017-10

**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 5 min

**Carrier Gas:** Constant Flow Helium, 0.8 mL/min

**Injection:** On-Column :1 0.5 μL @ 43°C

**Detection:** Mass Selective (MSD) (335°C)


**Analyst Note:** Analytes 10 ppm in Acetone



Products used in this application:



17589

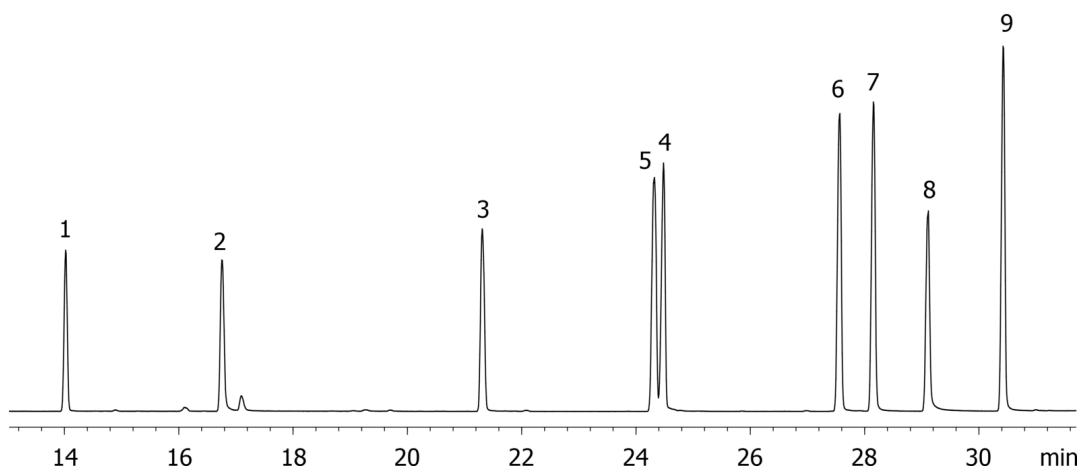

### ANALYTES:

- 1 EPTC
- 2 Mevinphos
- 3 Ethoprop
- 4 Atrazine
- 5 Propazine
- 6 Prometryn
- 7 Terbutryn
- 8 Triadimefon
- 9 Diphenamide

