


Opiates by MSD on a ZB-1MS**Column:** ZebtronTM ZB-1MS, GC Cap. Column 10 m x 0.18 mm x 0.18 μ m, Ea**Phase:** 100% Dimethylpolysiloxane**Dimensions:** 10 meters x 0.18 mm x 0.18 μ m**Order No:** 7CD-G011-08**Oven Profile:** 190 °C to 340 °C @ 20 °C/min**Carrier Gas:** Constant Flow Helium, 55 cm/sec**Injection:** Split 10:1 0.1 μ L @ 280°C**Detection:** Mass Selective (MSD) (230°C)**Analyst Note:** Analytes are 24ppm derivatized with BSTFA
Sim ions:
(Quant ion*) D6-Codeine (*377,349): Codeine (*371,343,234): D3-Morphine (*432,417): Morphine (*429,414,401)**Products used in this application:****ANALYTES:**

- 1 Codeine-d6
- 2 Codeine
- 3 Morphine-d3
- 4 Morphine

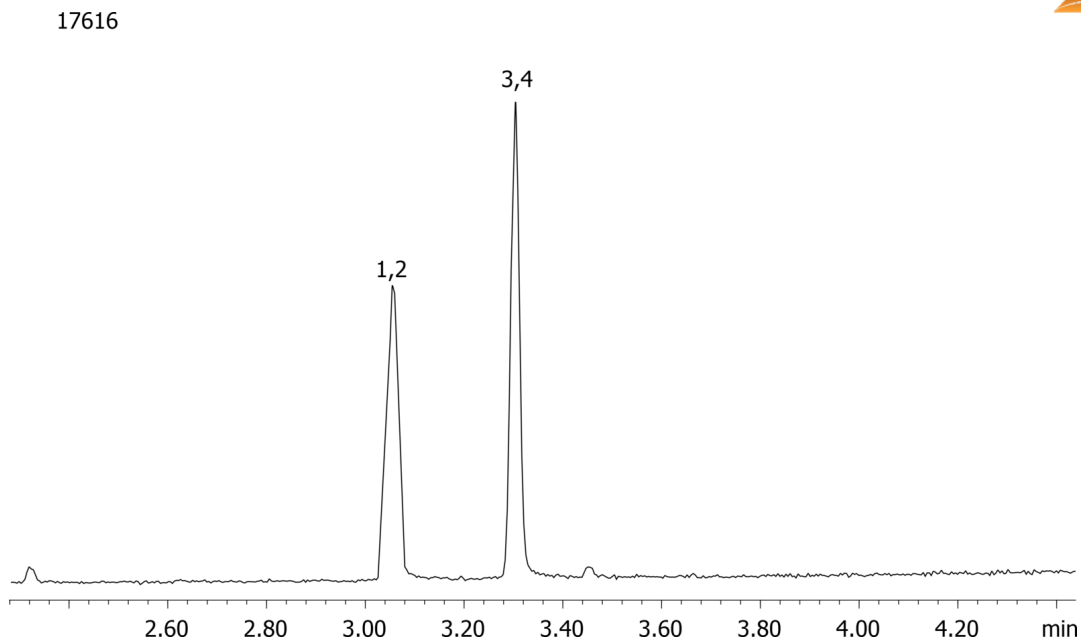