

# GC Application


ID No.: 17812

## Organophosphorous Pesticides Mix 2 on ZB-50

**Column:** Zebtron™ ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea  
**Phase:** 50% Phenyl 50% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25 µm  
**Order No:** 7HG-G004-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min  
**Carrier Gas:** Constant Flow Helium, 0.7 mL/min  
**Injection:** On-Column :1 0.5 µL @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes are 10 ppm in Acetone



Products used in this application:

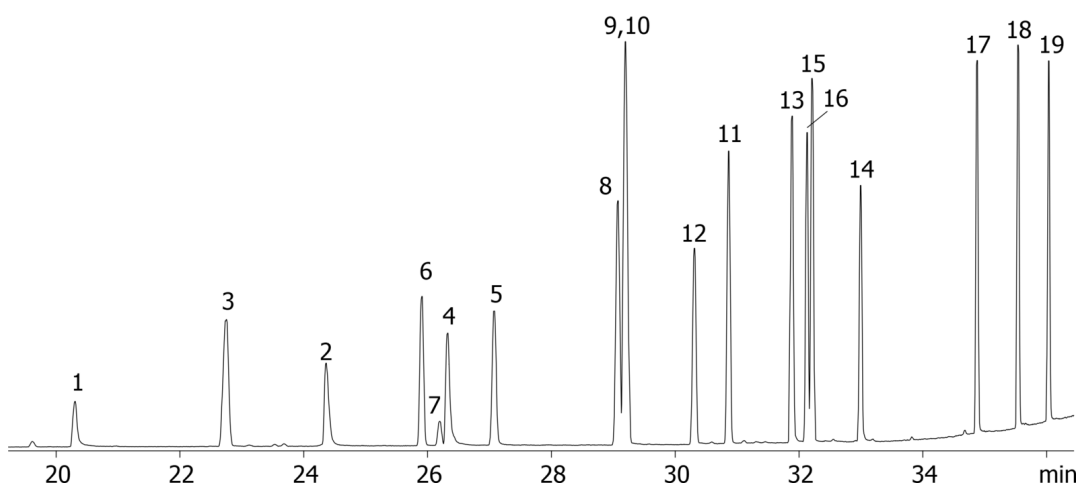

## Organophosphorous Pesticides Mix 2 on ZB-50

### ANALYTES:

- 1 Acephate
- 2 Omethoate
- 3 Cadusafos
- 4 Monocrotophos
- 5 Dimethoate
- 6 Tris(chloropropyl) phosphate
- 7 Tris(chloropropyl) phosphate
- 8 Tolclofos-methyl
- 9 Methyl parathion
- 10 Pirimiphos-methyl
- 11 Fenthion
- 12 Parathion
- 13 Quinalphos
- 14 Methidathion
- 15 Prothiophos
- 16 Tribufos
- 17 Triazophos
- 18 EPN
- 19 Phosalone

