


**Halogenated Pesticides Mix # 2 on ZB-50**

**Column:** Zebron<sup>™</sup> ZB-50, GC Cap. Column 30 m x 0.25 mm x 0.25  $\mu$ m, Ea  
**Phase:** 50% Phenyl 50% Dimethylpolysiloxane  
**Dimensions:** 30 meters x 0.25 mm x 0.25  $\mu$ m  
**Order No:** 7HG-G004-11  
**Oven Profile:** 40 °C for 1 min to 100 °C @ 15 °C/min to 225 °C @ 5 °C/min to 320 °C @ 15 °C/min for 8 min  
**Carrier Gas:** Constant Flow Helium, 0.7 mL/min  
**Injection:** On-Column :1 0.5  $\mu$ L @ 43°C  
**Detection:** Mass Selective (MSD) (335°C)  
**Analyst Note:** Analytes are 10 ppm in Acetone

**Products used in this application:**

17813

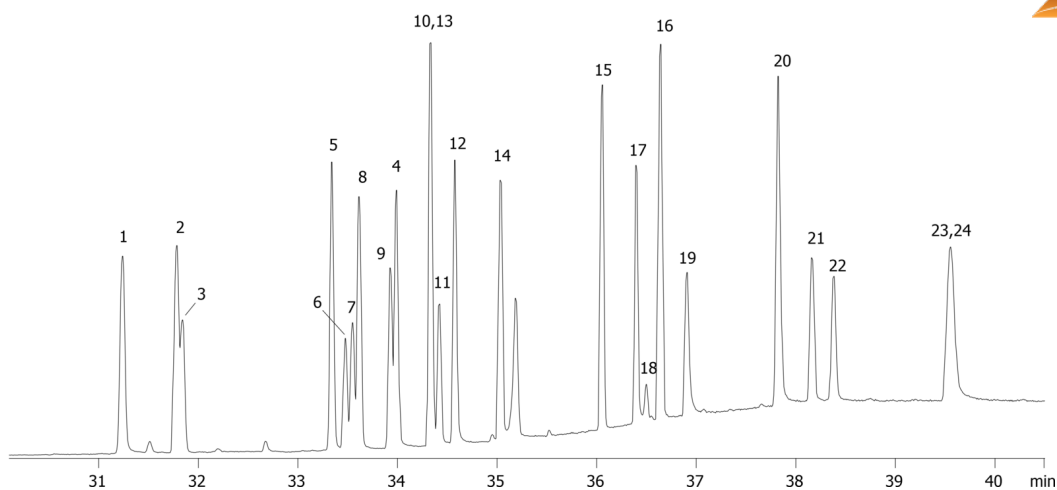

## Halogenated Pesticides Mix # 2 on ZB-50

### ANALYTES:

- 1 Penconazole
- 2 Procymidone
- 3 Endosulfan I
- 4 Fludioxonil
- 5 Myclobutanil
- 6 Cyproconazole
- 7 Cyproconazole
- 8 o,p'-DDT
- 9 Endosulfan II
- 10 Propiconazole
- 11 Propiconazole
- 12 Tebuconazole
- 13 Bifenthrin
- 14 Bromopropylate
- 15 Tetradifon
- 16 Fenarimol
- 17 cis-Permethrin
- 18 trans-Permethrin
- 19 Prochloraz
- 20 Fenbuconazole
- 21 Fenvalerate
- 22 Fenvalerate
- 23 Difenoconazole
- 24 Difenoconazole

